


Canon

Full Auto Tonometer

TX-20

Operation Manual

Make sure you read this manual before using the instrument.
Keep this manual in a safe place so that you can use it in the future.

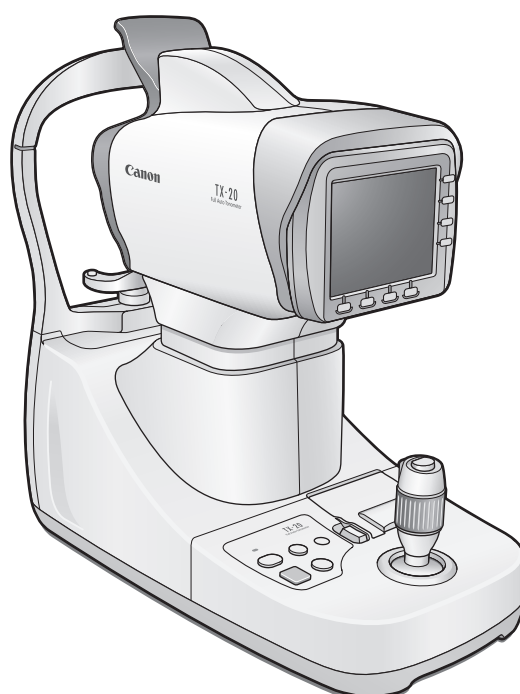

CE 0197

Thank you for purchasing the Canon Full Auto Tonometer TX-20 (referred to as “TX-20” in this manual). Be sure to read this manual thoroughly before using the instrument, and apply the information that you learn.

Important

- The TX-20 must only be used by a doctor or a legally qualified person.
- The user is responsible for managing the usage and maintenance of medical equipment. We suggest that a dedicated individual is assigned responsibility for maintenance to ensure that the TX-20 is kept in good condition and can be used safely.
- Rx Only (USA) Federal law restricts this device to sale by or on the order of a physician.

Disclaimers

- Canon takes no responsibility for damage that occurs due to fires, earthquakes, third-party actions, other accidents, the user's deliberate misuse, negligence, experimentation, or use under other abnormal conditions.
- Canon takes no responsibility for direct or consequential damages resulting from the use or the inability to use the TX-20.
- Canon takes no responsibility for injuries or property damage that may occur if safety precautions are not followed or the instrument is used for something other than its intended purpose.
- Medical examinations are the responsibility of a doctor. Canon takes no responsibility for diagnostic results.
- The information in this manual may be changed without prior notice.
- Although we have made every effort to ensure the accuracy of the information in this manual, if you have any questions regarding the contents, please contact your sales representative or local Canon dealer.

Installation

- Request your sales representative or local Canon dealer to install the TX-20.

Trademarks

- Canon and the Canon logo are trademarks of Canon Inc.
- Other system names and product names that appear in this manual are trademarks of their respective companies.

Copyright

- The copyright of this manual belongs to Canon Inc.
- Unauthorized reproduction, duplication, or reprinting of this manual in whole or in part is prohibited.

Contents

1 Introduction.....	5
Features.....	5
Indications for Use	5
Components.....	6
Conventions Used in This Manual	7
2 Safety	8
Regulatory Information.....	8
Safety Precautions	10
Notes on Using TX-20	14
Labels and Symbols.....	17
3 Names of Parts	19
Main Unit (From the Examiner's Side).....	19
Main Unit (From the Patient's Side).....	20
Bottom Side.....	21
Operation Panel.....	22
Operation Lever.....	23
Monitor Panel	24
Measurement Screen.....	25
PATIENT ID Screen.....	27
DATA Screen.....	28
MENU Screen.....	29
Screen that Appears When Measurement Fails	30
Anterior Segment Image Display Screen	30
Anterior Segment Image and Intraocular Pressure Measurement	
Waveform Display Screen.....	31
Other Screens	32
USB Connection Check Screen	32
Print Message Input Screen	32
User Name Input Screen	33
Password Input Screen	33
Computer Name Input Screen.....	34
Output Folder Name Input Screen	34
LCD Brightness Adjustment Screen	35
Description of Printouts	36
4 Basic Operation	37
Flow of Operations	37
Connecting Cables	38
Preparations Prior to an Examination	40
Patient Setup.....	41
Measuring.....	43
Finishing an Examination	47
Power Saving Mode	48
Recovering from the Power Saving Mode	48
5 Types of Measurements.....	49
Measurement Screen.....	49
Changing the Alignment Mode	49
Changing the Number of Measurements	50
Changing the Eye Fixation Lamp	51
PATIENT ID Screen	51
Input Patient ID in the PATIENT ID Screen	52
Input Patient ID from an External Device	56
MENU Screen	57
MENU Screen Display and Functions of Buttons	57

Basic Operations on MENU Screen	58
Common Settings in Measurement Modes	61
Settings for Intraocular Pressure Measurements	62
Settings for Printing	63
Settings for Patient ID Input Interface	64
Settings for Interface of Output of Measurements	64
Settings for Patient No.	67
Settings for Patient ID	68
Setting the Time, Date, and Power Saving Mode.....	69
Settings for the LCD	69
Checking Firmware	70
6 Maintenance.....	71
Daily Inspections	71
Cleaning and Disinfection	72
Refilling Chin Rest Paper	74
Replacing Printing Paper	74
Preparations for Packing and Transport	76
7 Troubleshooting.....	77
Symptom and Remedy	77
Message and Remedy	78
Appendix.....	87
Specifications.....	87
Stylesheet File.....	88
EMC (Electromagnetic Compatibility).....	90
Warranty and Repair Service	95

1 Introduction

Features

The Canon Full Auto Tonometer TX-20 has full auto-alignment for intraocular pressure measurement. The features of these products are shown below.

Three alignment modes

- In full auto mode, both eyes are measured and the results are printed automatically, shortening the time used for measuring.
- In auto mode, auto-alignment and measurement are performed each time the measurement button is pressed. The examiner can monitor the condition of the patient's eyes while measuring them.
- In manual mode, measurements can be performed manually such as when the patient's pupil is eccentric or cornea is deformed.

Compact, with superior operability

More compact and more lightweight than the previous model. Upright operations, patient care, alignment operation have been greatly improved with a tiltable color LCD monitor and powered operation lever.

Enhanced external interface

Measurement data can be exported via RS-232C connections or to a LAN. Patient ID can be input from a USB device.

Indications for Use




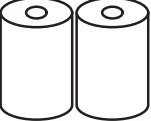



For USA

The Canon Full Auto Tonometer TX-20 is intended to be used for the measurement of intraocular pressure of the human eye.

For European Union

The TX-20 measures the intraocular pressure detecting the change of the cornea surface and the result will be provided for diagnosis.



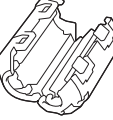



Components

Full Auto Tonometer TX-20 Main unit		Chin rest paper 100 sheets provided.	
Power cord Connects the TX-20 to an AC outlet (3m). USA and Canada only: Plug type: VM0275 Hospital-Grade		Printing paper 2 rolls provided.	
Nozzle cap Initially covers the nozzle.		Dust cover Use it to cover the TX-20 when the instrument is not being used.	
Ferrite core If a LAN cable is being used, attach the ferrite core to it.		Operation manual for the TX-20—this document Describes the handling precautions and operating instructions for the TX-20. Warranty Card (for USA model only) WEEE Directive Leaflet (for EU model only) Installation Report	

Optional products
Chin rest paper—500 sheets
Printing paper—10 rolls

Conventions Used in This Manual

This manual uses the following symbols to indicate precautions that are important for using the TX-20 safely. Always follow these safety precautions.

 WARNING	A warning that incorrect operation may result in death or serious injury.
 CAUTION	A caution that incorrect operation may result in serious injury.
CAUTION	A caution that incorrect operation may break the TX-20 or damage other devices.
	This symbol indicates actions that must not be taken (prohibited actions).
	This symbol indicates actions that must be taken.
	This symbol indicates important advice that we strongly recommend be followed while operating the TX-20.
	This symbol indicates supplementary explanations or advice for operating the TX-20.

References

This manual uses the following style to indicate reference destination.

Example: (see page 15)

2 Safety

Regulatory Information

The following sections list the classifications applicable to the TX-20 and the directives and standards that the TX-20 complies with.

Device Classification

Type of protection against electric shock	Class I equipment
Degree of protection against electric shock	Type B applied parts (Chin rest, forehead rest)
Degree of protection against ingress of water	IPX0
Degree of safety of application in the presence of a flammable anaesthetic mixture with air or with oxygen or nitrous oxide	Not suitable
Mode of operation	Continuous operation

Directives and Standards

USA and Canada

UL 60601-1:2003 CAN/CSA C22.2 No.601.1-M90 IEC 60601-1: 1988/A1: 1991/A2: 1995	Medical electrical equipment – Part 1: General requirements for safety
IEC 60601-1:2005	Medical electrical equipment - Part1: General requirements for basic safety and essential performance
IEC 60601-1-1: 2000	Medical electrical equipment – Part 1-1: General requirements for safety – Collateral standard: Safety requirements for medical electrical systems
IEC 60601-1-2: 2001/A1:2004 IEC 60601-1-2:2007	Medical electrical equipment – Part 1-2: Collateral standard: Electromagnetic compatibility – Requirements and tests
IEC 60601-1-4: 1996/A1: 1999	Medical electrical equipment – Part 1-4: Collateral standard: Programmable electrical medical systems
IEC 62304: 2006	Medical device software – Software life-cycle processes
IEC 62366: 2007	Medical device – Application of usability engineering to medical devices
ISO 15004-1: 2006	Ophthalmic instruments – Part 1: General requirements applicable to all ophthalmic instruments
ISO 15004-2: 2007	Ophthalmic instruments – Part 2: Light hazard protection

ISO 10993-1: 2009
ISO 10993-5: 2009
ISO 10993-10: 2010

Biological evaluation of medical devices – Evaluation and testing

ISO 8612: 2001
ANSI Z80.10-2003

Ophthalmic instruments – Tonometers

European Union

93/42/EEC

Medical Device Directive

EN 60601-1:2006

Medical electrical equipment - Part1:General requirements for basic safety and essential performance

EN 60601-1-2: 2007

Medical electrical equipment – Part 1-2: Collateral standard: Electromagnetic compatibility – Requirements and tests

EN 60601-1-6: 2010

Medical electrical equipment – Part 1-6: Collateral standard: Usability

EN 62304: 2006

Medical device software – Software life-cycle processes

EN 62366: 2008

Medical device – Application of usability engineering to medical devices

EN ISO 14971: 2009

Medical device – Application of risk management to medical devices

EN ISO 15004-1: 2009

Ophthalmic instruments – Part 1: General requirements applicable to all ophthalmic instruments

ISO 15004-2: 2007

Ophthalmic instruments – Part 2: Light hazard protection

EN ISO 10993-1: 2009
EN ISO 10993-5: 2009
EN ISO 10993-10: 2010

Biological evaluation of medical devices

ISO 8612: 2001

Ophthalmic instruments – Tonometers

Safety Precautions

To prevent injuries and data loss, operate the TX-20 correctly by following the safety precautions.

WARNING



To avoid the risk of electric shock, this instrument must only be connected into a supply mains with protective earth.



Prohibited

Do not damage the power cord.

- Do not place anything heavy onto the power cord.
- Do not damage or modify the power cord.
- Do not forcibly bend, twist or pull the power cord.
- Do not hold the power cord when removing it from the AC outlet. Be sure to hold the plug.

Handle the power cord carefully. If the cord is damaged, contact your sales representative or local Canon dealer for a replacement. A damaged cord may result in fire or electric shock.



Prohibited

Do not use a multiple power strip and an extension cable.

Connect the power cord directly to the AC outlet. Do not use a multiple power strip and an extension cable with it.



Prohibited

Do not disassemble or modify.

A disassembled or modified instrument may result in fire or electric shock. Because the TX-20 incorporates high-voltage parts that may cause electric shocks, touching them may cause death or serious injury.



Prohibited

Do not leave alcohol, thinner, or any flammable chemicals near the instrument.

Do not place near to a flammable solvent. Fire may result if the solvent spills or evaporates and makes contact with the internal electric parts. Some disinfectants are flammable. Be sure to take precautions when using them.



Stop using immediately if there is an abnormality or problem. If an abnormality occurs, immediately unplug the power plug and turn off the power of all the devices.

- Smoke is emitted
- An unusual smell
- An unusual noise
- Foreign matter gets inside
- Damaged

Fire or electric shock may result if you continue use in such conditions. Immediately turn off the power of TX-20, unplug the power plug and turn off the power of all connected devices. Then, contact your sales representative or local Canon dealer.



Prohibited

Do not place anything on top of the device.

Fire or electric shock may result if water or any other liquid or a needle, paper clip or any foreign matter gets inside the TX-20.

⚠ WARNING

Prohibited

Do not use a power supply voltage other than the one specified on the rating label.**Do not use any power cord other than the one supplied.**

Use the power voltage specified on the rating label. Otherwise, fire or electric shock may result.

The supplied power cord is designed exclusively for this product. Do not use any other power cord.



Prohibited

Do not plug or unplug the power plug with wet hands.

Do not plug or unplug the power plug or handle any other parts with wet hands. Otherwise, fire or electric shock may result.

**Unplug the plug periodically and remove any dust or dirt around the plug and the AC outlet.**

If the cord is kept plugged in for a long time in a dusty, humid or sooty environment, dust around the plug will attract moisture, resulting in possible insulation failure that may result in a fire.

**Insert the power plug completely.**

Insert the power plug completely into the AC outlet. If a pin of the power plug makes contact with metal or any conductive object, fire or electric shock may result.



Prohibited

Do not clean the TX-20 with flammable solvent.

For safety reasons, before cleaning the TX-20, be sure to turn off the power of all the devices and unplug the power cord from the AC outlet. Do not use alcohol that is not for cleaning, benzene, thinner or any other flammable solvent. Otherwise, fire or electric shock may result.

**Turn off the power before inspection.**

For safety reasons, before inspecting the instrument or cables, be sure to turn off the power of all the devices. Otherwise, electric shock may result.



Prohibited

Do not measure an eye whose cornea is frail due to corneal disease or surgery.

Otherwise, complications may occur.

**Be sure to turn off the power before moving the TX-20.**

When moving the TX-20, always move the measurement unit to the packing position (see page 76), turn off the power, unplug the power cord, and disconnect any external equipment. Otherwise, the cable may be damaged, which may result in fire or electric shock.

**Use the indentations for lifting and moving the TX-20.**

When lifting the TX-20, hold it on the bottom and keep it level. Do not hold it by the head rest, nozzle or LCD monitor, as they may come off and result in injury.



Prohibited

Do not touch conductive parts of non-medical equipment and the patient simultaneously.

Otherwise, electric shock may result.

⚠ CAUTION

Prohibited

Do not install in locations exposed to water, steam, moisture or dust.

Doing so may cause problems or malfunctions.

CAUTION**Do not install in locations exposed to salt, sulfur or corrosive gas.**

Doing so may result in corrosion of the instrument, problems or malfunctions.

**Do not install in locations that are unstable or exposed to vibration.**

The vibration may knock over the instrument or the instrument may become unbalanced and fall, resulting in a malfunction or injury.

**Do not place anything near the power plug.**

To make it easy to unplug the power plug at any time, avoid putting any obstructions near the AC outlet.

Fire or electric shock may result if the power plug is not unplugged during an emergency.

**Do not block the vent holes.**

Doing so may cause the internal temperature to rise and result in fire.

**Do not put your hands or fingers near the measurement unit.**

Your hand or fingers may be pinched and injured.

Similarly, instruct the patient not to place his/her hands or fingers on the base.

**Hold the TX-20 when connecting or disconnecting a cable.**

For safety reasons, when connecting or disconnecting the power cord or any cable be sure to hold the main unit. Otherwise, the main unit may fall over, possibly causing injury.

**Ensure that the entire system conforms to IEC 60601-1-1.**

Use the computer and other equipment that conform to the system standard IEC 60601-1 or IEC 60950-1 for the Full Auto Tonometer TX-20. Be sure that the entire system conforms to IEC 60601-1-1. Be sure to also use an isolation transformer conforming to IEC 60601-1 when a computer conforming to IEC 60950-1 is used. Otherwise, electric shock may occur. Ask your sales representative or local Canon dealer to connect it to a computer.

**Use an isolation transformer for LAN connections.**

If the LAN connector is being used to connect to a network, connect an isolation transformer for networks between the TX-20 and network devices (HUB etc.) and between network devices.

If this is not done, there is a risk of short circuits or malfunctioning devices. Ask your sales representative or local Canon dealer to connect it to a network.

**Turn off the power and unplug the power cord before connecting cables.**

Always turn off the power and unplug the power cord before connecting the cables of any external devices to the TX-20.

If this is not done, there is a risk of damage.

**Always disinfect the nozzle after measuring an infected patient's eyes.**

A patient's tears may get on the nozzle and cause a secondary infection. Always use a cotton swab or other appropriate material that has been treated with disinfectant ethanol to wipe the nozzle.

CAUTION**Keep the forehead rest and chin rest clean.**

To prevent the risk of infection, wipe the forehead rest with disinfectant ethanol for each patient. To ensure cleanliness, replace the chin rest paper for each patient. If the chin rest paper is not being used, be sure to disinfect the chin rest for each patient.

For details on how to disinfect, consult a specialist. The forehead rest may be corroded if a disinfectant other than those above is used.

**Be sure to turn off the power when not in use.**

For safety reasons, be sure to turn off the power of all the devices when the TX-20 is not being used.

Also, remove the power plug from the AC outlet and put on the cover when the TX-20 is not going to be used for a long time.

Otherwise, dust or any foreign matter may accumulate and result in fire.

**Make sure to select the position where the nozzle stops for each patient.**

Look from the side of the patient when setting the stop position to see that the patient's forehead and chin are firmly against the forehead rest and chin rest. If the stop position is not set correctly, the nozzle may contact the patient's eye and result in injury.

**Prohibited****Do not touch the cutter for the printer.**

Do not touch the cutter for the printer. Similarly, instruct the patient not to touch the cutter for the printer.

Otherwise, it may result in injury.

**Prohibited****Do not touch the area around the printer or open the paper cover during or immediately after printing.**

The printer's thermal head and the surrounding area are very hot during and immediately after printing. To prevent burns, wait for the printer to cool before opening the paper cover to replace paper.

**Prohibited****Do not touch the measurement unit and the chin rest while they are moving.**

The measurement unit and the chin rest move when the power is turned on, full auto is set, or the CHIN REST button is pressed. Do not touch the measurement unit or the chin rest while they are moving. Keep the patient's chin away from the chin rest.

CAUTION**Use the product's packaging to transport it.**

Use the TX-20 packing materials to protect it from vibration and shocks when it is transported on a cart, in a vehicle, or if shipped over long distances.

Vibration or shock may cause a breakdown or damage the TX-20.

For details, contact your sales representative or local Canon dealer.

**Inspect daily and periodically.**

For safety reasons, before using the TX-20, be sure to perform the daily inspection.

Have a periodically inspection performed for the TX-20 at least once a year by a Canon designated representative to maintain its performance and reliability.

Notes on Using TX-20

Before Use

- Sudden heating of the room in cold areas will cause condensation to form on the optical component of the nozzle and on optical parts inside the instrument. In this case, wait until condensation disappears before performing a measurement.

During Measurement

- Describe to the patient what to expect during measurement so he/she will not be surprised by the sudden air puff. Demonstrate the intensity of the air puff by having the patient place a finger in front of the nozzle and then press the CHECK button.
- The measurement unit automatically moves after the power is turned on. Do not prevent movement of the measurement unit by holding it or placing something near it. Do not let the patient place his/her chin on the chin rest until the measurement unit stops.
- It is recommended that a hard copy of the printout be made if you wish to store it for a long time, because printouts on thermal paper are apt to deteriorate.

After Use

- Turn off the power, put the cap over the optical component of the nozzle in order to prevent dust from attaching to it, and place the dust cover over the instrument.

During Cleaning and Disinfection

- If the surface of the optical component of the nozzle is wiped when dust or dirt is on it, it will be scratched. Also, do not wipe the optical component with alcohol, including disinfectant ethanol or silicone-coated paper. Otherwise, the surface will be corroded or stained.
- Use neutral detergent for cleaning the outside of the TX-20. Do not use neutral detergent for cleaning the LCD monitor. See below for LCD cleaning instructions.

LCD Monitor

- Approximately 99.99% of the dots in high-technology LCD monitors are functional. About 0.01% of the dots may be missing or might appear as a bright spot. This is not a malfunction.
- The polarization plate on the surface of the monitor is very fragile. Be careful when handling it and do not press or rub any hard objects against it.
- Wipe water or saliva from the polarization plate on the monitor's surface immediately; leaving it for a long time may cause distortion or discoloration.
- If the surface is dirty, use something soft, such as an absorbent cotton cloth, with a small amount of ethanol, and lightly wipe it.

Environment of Use

- Use, store, and transport the TX-20 in an environment that is within the range of the following conditions. Use the original packaging to store or ship it.

	Temperature	Humidity	Atmospheric pressure
Environment of use	10 to 40°C	30 to 90% RH (no condensation)	600 to 1060 hPa
Storage and transportation environment	–30 to 50°C	10 to 95% RH (no condensation)	600 to 1060 hPa

- Do not install, store, or leave the TX-20 in a very hot or humid environment. Also, do not use the TX-20 outside. Doing so may cause a malfunction or misoperation.
- After many years of usage, airborne dust in the room may get on the nozzle as well as the optical parts in the main unit and prevent accurate measurement.
- When the TX-20 is not being used, attach the nozzle cap and place the dust cover over the TX-20.

Installation

- Request your sales representative or local Canon dealer to install the TX-20.
- A strong shock to the TX-20 may put it out of alignment. Please handle the unit carefully.

Transportation

- Be sure to turn off the TX-20's power switch, remove the AC plug from the outlet, and disconnect any cables connected to other equipment.
- When moving the TX-20, perform the specified operation to move the measurement unit to the transport position (see page 76).
- Remove the paper from the printer before moving the TX-20.
- Carry the TX-20 by the indentations for lifting at its front and rear, and keep it level.
- Do not hold anything except the indentions for lifting when lifting the TX-20.
- The TX-20 needs to be protected from vibration and shocks when it is transported on a cart, in a vehicle, or if shipped over long distances. Always use the original packaging. For details, contact your sales representative or local Canon dealer.

Disposal

Disposal of this product in an unlawful manner may have a negative impact on human health and on the environment. Therefore, when disposing of this product, be absolutely certain to follow the procedure which is in conformity with the laws and regulations applicable to your area.

European Union (and EEA*) only



This symbol indicates that this product is not to be disposed of with your household waste, according to the WEEE Directive (2002/96/EC) and your national law. This product should be handed over to a designated collection point, e.g., on an authorized one-for-one basis when you buy a new similar product or to an authorized collection site for recycling waste electrical and electronic equipment (EEE). Improper handling of this type of waste could have a possible negative impact on the environment and human health due to potentially hazardous substances that are generally associated with EEE. At the same time, your cooperation in the correct disposal of this product will contribute to the effective usage of natural resources. For more information about where you can drop off your waste equipment for recycling, please contact your local city office, waste authority, approved WEEE scheme or your household waste disposal service. For more information regarding return and recycling of WEEE products, please visit www.canon-europe.com/WEEE.


* EEA: Norway, Iceland and Liechtenstein

For California, USA Only

The battery in this device contains perchlorate material. Special handling may apply. See <http://www.dtsc.ca.gov/hazardouswaste/perchlorate>

Labels and Symbols

The position and contents of the label attached on the TX-20 are shown below.

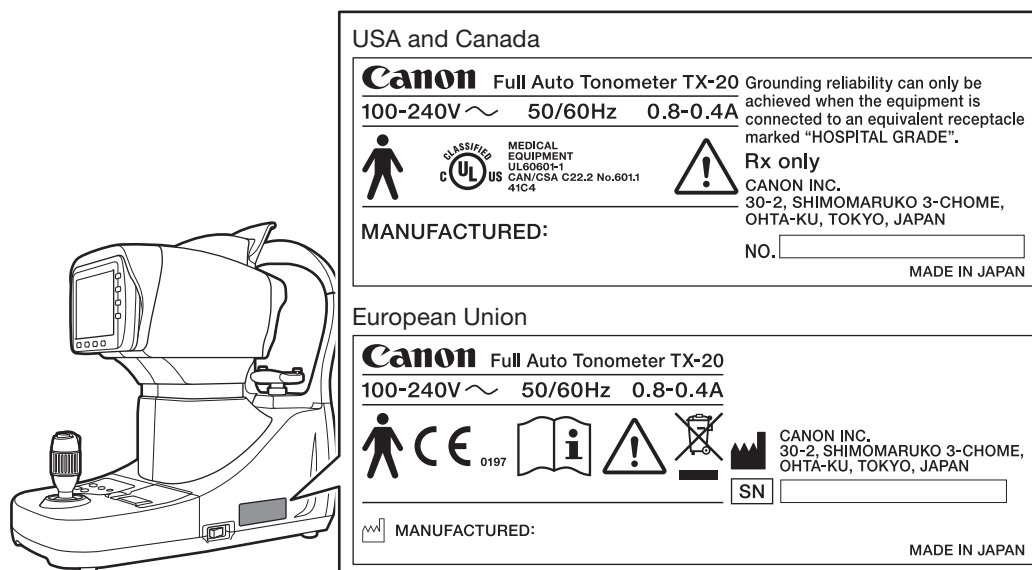










The following table describes the symbols and indications on the TX-20.

USA and Canada


	Alternating current
	Type B
	Caution: Check the documentation provided.
	MEDICAL EQUIPMENT WITH RESPECT TO ELECTRICAL SHOCK, FIRE AND MECHANICAL HAZARDS ONLY IN ACCORDANCE WITH UL 60601-1 AND CAN/CSA C22.2 NO. 601.1 41C4
Rx only	Caution: Federal law (USA) restricts this device to sale by or on the order of a physician.
MANUFACTURED:	Year and month of manufacture Example: October 2010
NO.	Serial number in SIX digits Example: 123456

European Union

	Alternating current	
	Type B	
	Caution: Check the documentation provided.	
	Product that WEEE directive, Directive on Waste Electrical and Electronic Equipment, 2002/96/EC, requires separate collection.	
	This symbol shows compliance of the equipment with Directive 93/42/EEC.	
	This symbol indicates that a User's Manual is supplied with this equipment.	
	Manufacturer's name and address	
	Serial number in SIX digits	Example: 123456
	Year and month of manufacture	Example: October 2010


3 Names of Parts

Main Unit (From the Examiner's Side)

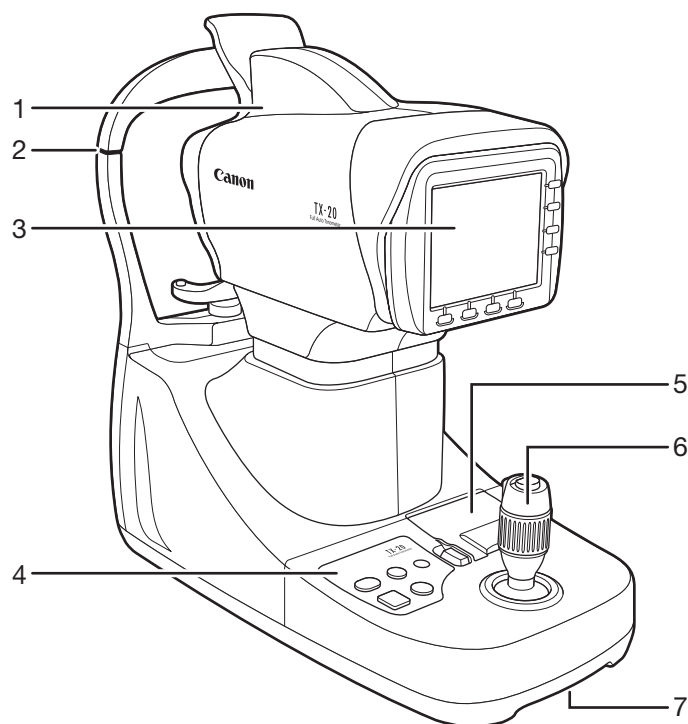


- | | | | |
|---|--|---|--|
| 1 | Measurement unit
Measurement is done at this part. | 5 | Printer
Prints the measurement results. |
| 2 | Height adjustment mark
Align the patient's eye with this mark. | 6 | Operation lever |
| 3 | LCD monitor
The patient's eye and various screens appear here. | 7 | Indentation for lifting |
| 4 | Operation panel
Operations such as setting the nozzle position,
chin rest height, and printing are done from here. | | |

Main Unit (From the Patient's Side)

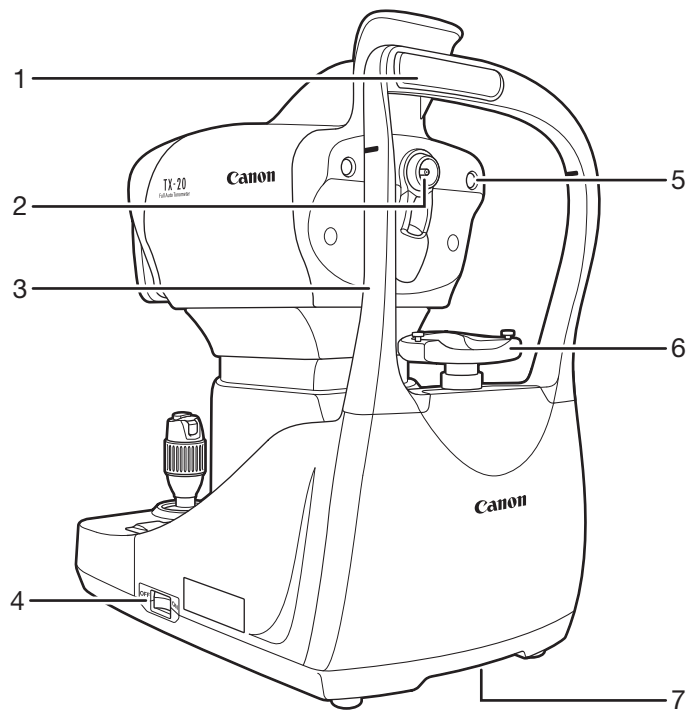

- | | | | |
|---|---|---|--|
| 1 | Forehead rest | 5 | External eye fixation lamp |
| | Patient rests his/her forehead here to prevent patient's eye from moving. | | Used if it is difficult to fix the patient's eye by using the internal eye fixation lamp. Fix on the eye that is not being measured. |
| 2 | Nozzle | 6 | Chin rest |
| | Blows air for observing or measuring the patient's eye. | | The patient's chin is placed here. Raise or lower the chin rest to set the appropriate height for the patient's eye(s). |
| 3 | Head rest | 7 | Indentation for lifting |
| 4 | Power switch | | |
| | Turns on the power for the main unit. | | |

Bottom Side

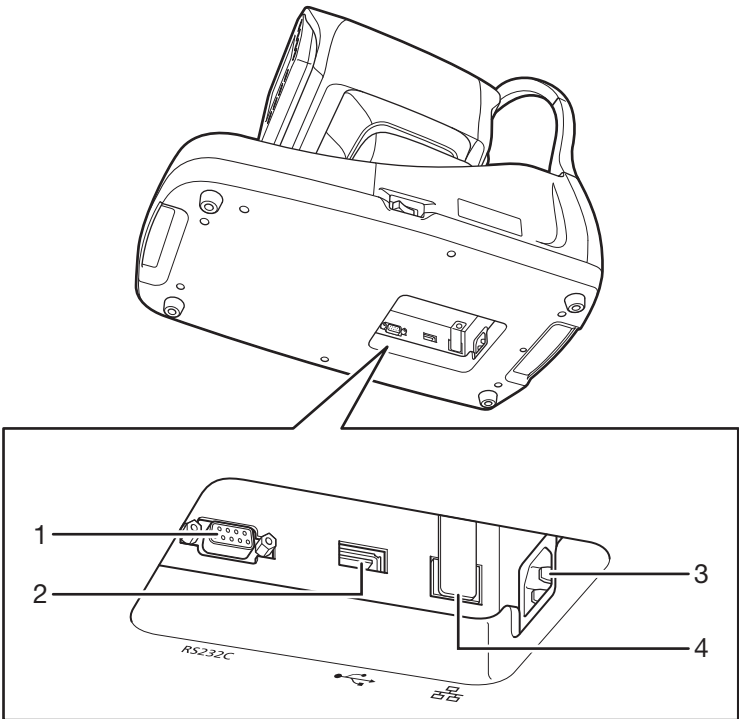


- | | | | |
|---|--|---|---|
| 1 | RS-232C connector | 3 | AC power connector |
| | Transfers measurement data to an external device. | 4 | LAN connector |
| | Can be connected to a device that supports the RS-232C standard. | | Transfers measurement data to an external device. |
| 2 | USB connector | | |
| | Reads patient ID from an external device. | | |

!

Using the RS-232C or LAN connector

See “Safety Precautions” (see page 12) and “Settings for Interface of Output of Measurements” (see page 64).


Operation Panel

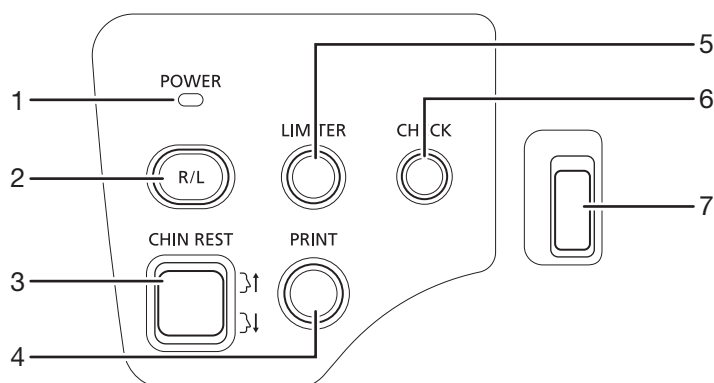

- | | |
|--|---|
| <p>1 POWER lamp
Lights when power is turned on.</p> <p>2 R/L button
Switches between the patient's left and right eye as the measurement unit moves automatically.</p> <p>3 CHIN REST button
Moves the chin rest up and down to align the height of the patient's eye.</p> <p>4 PRINT button</p> <ul style="list-style-type: none"> • Prints the measurement results. Sends the measurement results to an external device if one is connected and the settings have been done in advance. • Pressing the PRINT button for 2 seconds or more prints or transmits the measured value and standard value for the number of measurements (see page 63). • For the preparation for the next measurement, the settings change as shown below. <ul style="list-style-type: none"> - Nozzle stop position: Releases the setting. - Selection of eye fixation lamp: The internal eye fixation lamp is selected. <p>The external eye fixation lamp turns off and the internal eye fixation lamp flashes.</p> - Patient ID: The patient ID is cleared* (not entered state). - Patient No: Counts up the patient numbers (however, according to settings). - Display of measured value: Measured values are cleared* (not measured state). - Measurement unit: Returns to the initial position. | <p>5 LIMITER button
Sets the stop position to prevent the measurement unit from touching the patient.</p> <p>6 CHECK button
Blows air from the nozzle.</p> <p>7 Paper cover button</p> |
|--|---|

*You can confirm the patient ID or measured value on the DATA screen until the measurement button is pressed.

Operation Lever

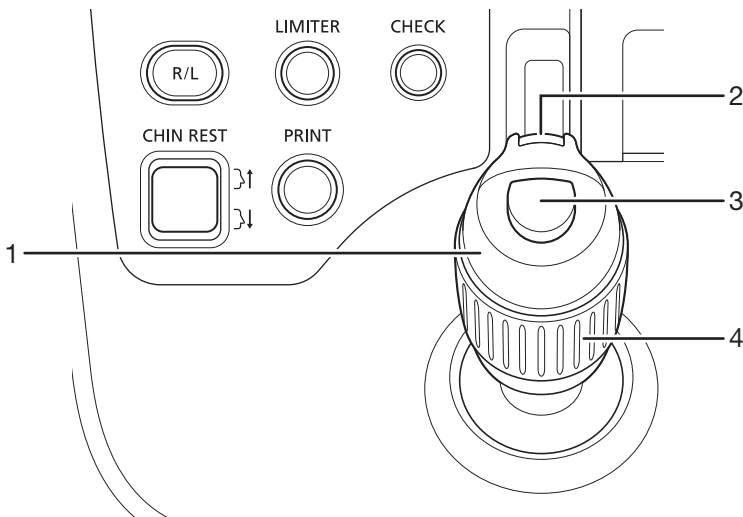

- 1

Operation lever

Tilt the operation lever to move the measurement unit.

Tilt right or left: Moves slightly to right or left.

Tilt front or back: Moves slightly to the front or back.

Tilt past the maximum to the front, back, left, or right: Motion is great.
- 2

Alignment off button

While the button is pressed, the measurement unit does not move even if the operation lever is tilted. Use to adjust the angle of the operation lever without moving the measurement unit. For example, it can be used if the operation lever is at an awkward and difficult-to-manipulate angle when aligning with the patient's eye.
- 3


Measurement button

Starts the measurement.
- 4


Measurement unit's vertical movement ring

Turn the ring right and left to move the measurement unit up and down.

Monitor Panel

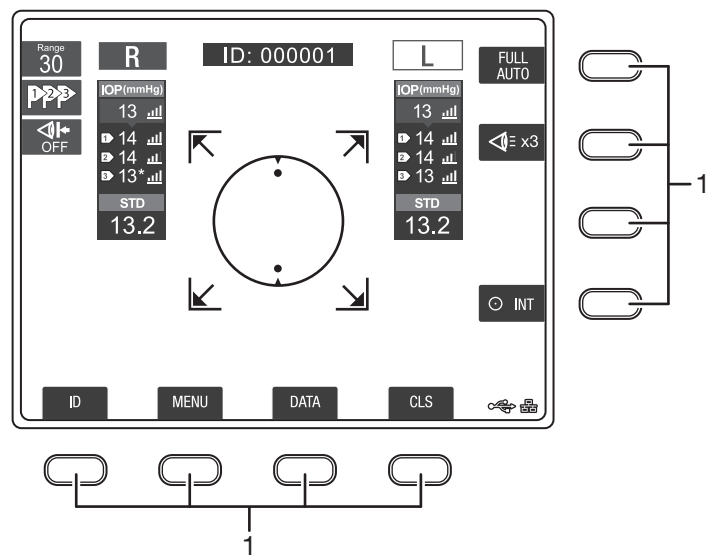


- 1 Soft keys
- The purpose of soft keys varies according to the assigned function. An icon on the monitor near the soft key indicates the purpose assigned to the soft key.
- Pressing a soft key changes and sets functions, and changes the screen.
- In this manual, these soft keys are indicated by the name of their functions.

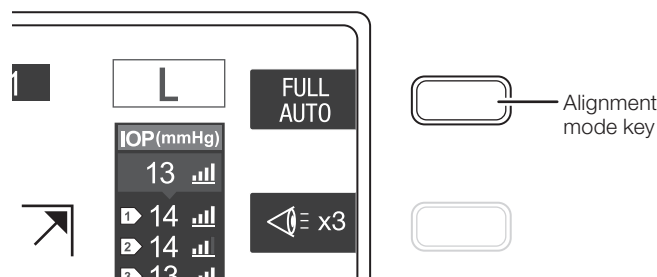

Measurement Screen

This screen appears when doing measurements.
You can set the measurement functions, show the most recent measurement data, and switch screens.

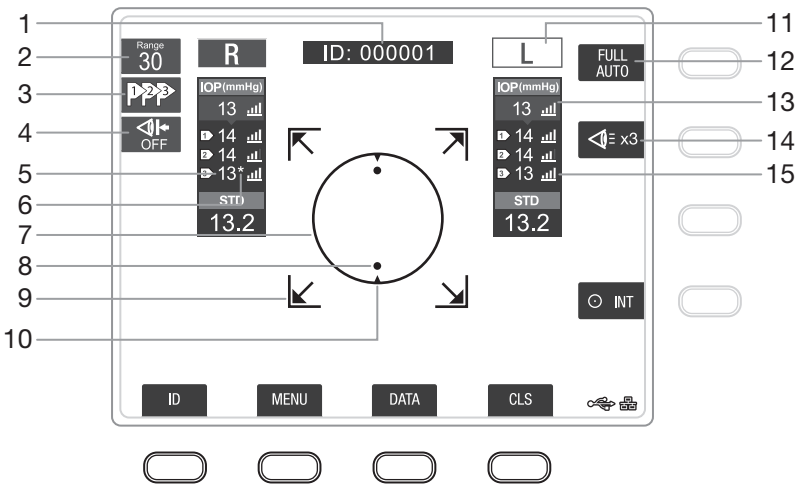

- 1

Patient ID/Patient No.
- 2

Measurement range
- 3

Display order of measurement value
1 2 3 : From highest to lowest reliability, 1 to 3.
1 2 3 : From newest to oldest measurement, 1 to 3.
- 4

LIMITER OFF mark
Indicates the LIMITER function is disabled.
- 5

Measurement results history
Shows the results of measurements according to age or reliability.
- 6

Low reliability mark
This mark appears when reliability of the measurement is considered low as the result of the analysis of the measured value.
- 7

Alignment circle
Mark for aligning vertical and horizontal position of patient's eye and measurement unit.
- 8

Corneal light point
A point that appears optically, indicating alignment of the position to measure in the patient's eye and the measurement unit.
- 9

Measurement unit motion limit mark
Indicates measurement unit has moved to the limit position.
- 10

Alignment guide
This guide is an index that shows the corneal light point. It only appears in manual alignment mode.
- 11

Left/right eye indicator
R L : Right eye
R L : Left eye
- 12

Alignment mode setting
Press the adjacent alignment mode key to switch the alignment mode.
FULL AUTO → AUTO → MANUAL
- 13


Most recent intraocular pressure
Shows the most recent measurement.
- 14

Measurement count setting
Press the adjacent measurement count key to switch the measurement count setting. Only enabled in full alignment mode.
x1 → x2 → x3
- 15

Level of received light
Indicates the amount of light received from light reflected from cornea when measuring intraocular pressure. The number of lines increases as the amount of received light increases.

If there is little received light, light reflected from the cornea may be obstructed by the eyelid or eyelashes. However, it is also affected by the corneal reflectivity of the patient.

3 Names of Parts

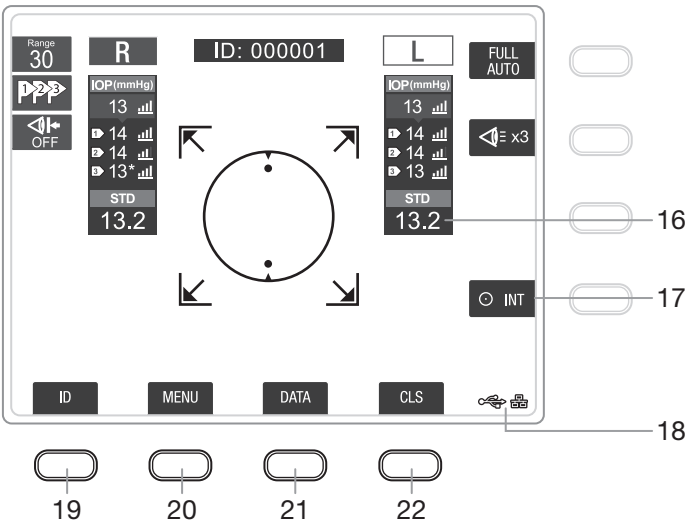

- 16

Intraocular pressure standard value

Representative intraocular pressure value that was statistically derived from measured values. Measured values with a low reliability mark are not used.
- 17

Eye fixation lamp

Press the adjacent eye fixation lamp key to switch the internal and external eye fixation lamps.

○ INT

→

○ EXT
- 18

External device connection status

Shows the status of external device connection.
- 19

ID key

Press to switch to the **PATIENT ID** screen.
- 20

MENU key

Press to switch to the **MENU** screen.
- 21

DATA key

Press to switch to the **DATA** screen.
- 22

CLS key


Press to do the following.

- Nozzle stop position: Releases the setting.
 - Selection of eye fixation lamp: The internal eye fixation lamp is selected.
The external eye fixation lamp turns off and the internal eye fixation lamp flashes.
 - Patient ID: The patient ID is cleared* (not entered state).
 - Patient No: Counts up the patient numbers (however, according to settings).
 - Display of measured value: Measured values are cleared* (not measured state).
 - Measurement unit: Returns to the initial position.

*You can confirm the patient ID or measured value on the DATA screen until the measurement button is pressed.

PATIENT ID Screen

This screen is used to input patient IDs.

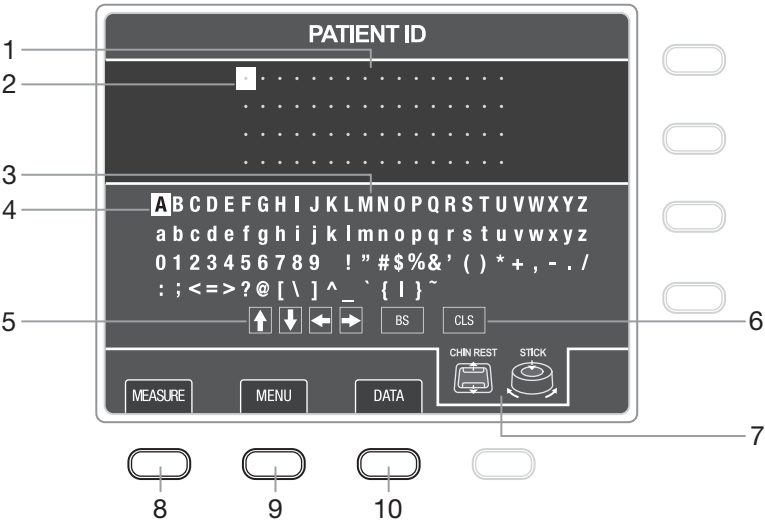

- | | |
|---|--|
| <p>1 Patient ID area
The text selected from the list of characters appears in the patient ID input field.
A maximum of 64 ASCII code characters can be input.</p> <p>2 ID cursor
Shows the position where the patient ID can be edited.</p> <p>3 List of characters
List of characters that can be input for a patient ID.</p> <p>4 Selection cursor
Select characters to input for a patient ID.</p> <p>5 ID cursor movement commands
The ID cursor moves in the direction of the selected arrow.</p> <p>6 ID delete command
[BS]: Deletes the character that is before the character specified by the ID cursor. However, if the ID cursor is specifying the first character, it deletes that character.

[CLS]: Deletes all the characters in the patient ID. The confirmation screen is displayed; select one of the following.</p> | <p>7 Operation reference guide
Shows the types of buttons that can be operated in the PATIENT ID screen. Refer to "5 Types of Measurements" for operating procedures.</p> <p>8 MEASURE key*
Press to switch to the Measurement screen.</p> <p>9 MENU key*
Press to switch to the MENU screen.</p> <p>10 DATA key*
Press to switch to the DATA screen.</p> |
|---|--|
- [OK]: Deletes all text.
 - [CANCEL]: Does not delete.

*By doing the settings in advance, if a patient ID has not been input when it is pressed, a message indicating that the patient ID has not been input is displayed (see page 68).

DATA Screen


Diagram illustrating the layout of the device screen and controls, with numbered callouts:

- 1: Screen area (R, DATA, L)
- 2: R button
- 3: L button
- 4: IOP (mmHg) header
- 5: IOP (mmHg) header
- 6: IOP (mmHg) header
- 7: ID button
- 8: MEASURE button
- 9: MENU button
- 10: DELETE button

- 28

MENU Screen

Perform the settings for measurements, printing, and operating the measurement unit.

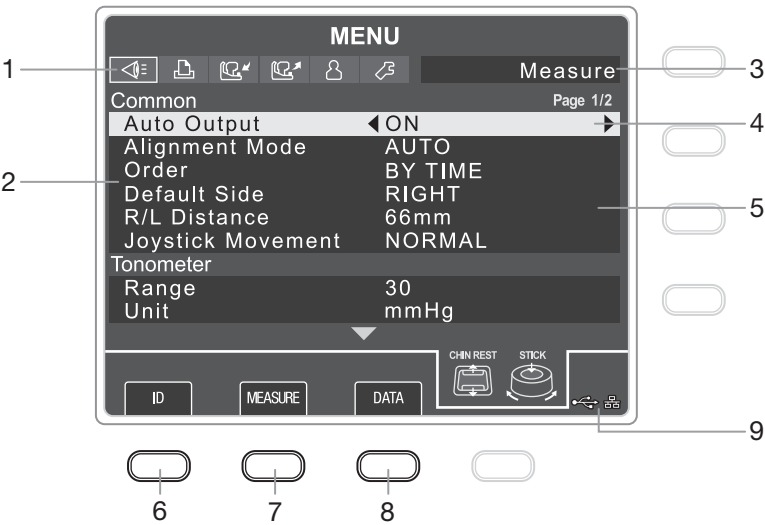

- 1

Settings tab

The items to set are grouped according to their purpose. Selected items are shown in orange.

Measure tab: Settings for intraocular pressure measurements

Print tab: Settings for printing with the built-in printer

Input tab: Settings for inputting data from an external device

Output tab: Settings for outputting data to an external device

Patient tab: Settings for patient ID and patient No.

Other tab: Settings for time, power saving mode, and the LCD, and shows information about the firmware.
- 2

Setting items

The setting items for each tab are displayed.
- 3

Name of the selected tab
- 4

Cursor

Shows the position of the selected setting value, setting item, and tab. They light orange when selected.
- 5

Settings

The settings for each setting item are displayed.
- 6

ID key

Press to switch to the **PATIENT ID** screen.
- 7

MEASURE key

Press to switch to the **Measurement** screen.
- 8

DATA key

Press to switch to the **DATA** screen.
- 9

External device connection status

Shows the status of external device connection.

: USB device is connected

: Error has occurred in USB device connection

: LAN is connected

: RS-232C device is connected


No display: No connection

Screen that Appears When Measurement Fails

Anterior Segment Image Display Screen

If a measurement error occurs while measuring intraocular pressure, the anterior segment image from just prior to measuring is displayed. The anterior segment image disappears after 5 seconds or if you press a soft key.

If [OFF] is set for [Error Display] in the intraocular pressure measurement menu, the anterior segment image is not displayed (see page 62).

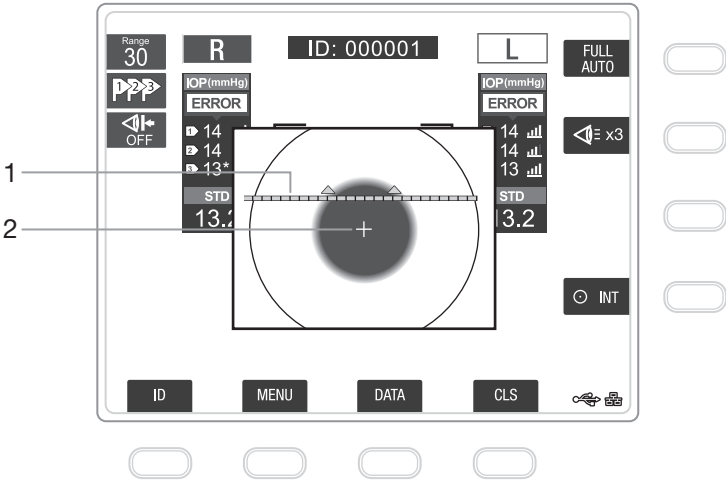

- 1

Eyelid warning line

Perform the measurement while the eyelid and eyelashes are above the warning line.
- 2

Cross mark


Displays the intraocular pressure measurement optical axis.

Anterior Segment Image and Intraocular Pressure Measurement Waveform Display Screen

If the values of intraocular pressure measurements are uniformly dispersed, the TX-20 displays the anterior segment images and the intraocular pressure waveforms that are captured just prior to each of the three latest measurements as below. This screen display is available for the measurements performed over three times.

Press the soft key to hide the anterior segment image and intraocular pressure measurement waveform.

If [OFF] is set for [Warning Display] in the intraocular pressure measurement menu, the anterior segment image and intraocular pressure measurement waveform are not displayed (see page 62).

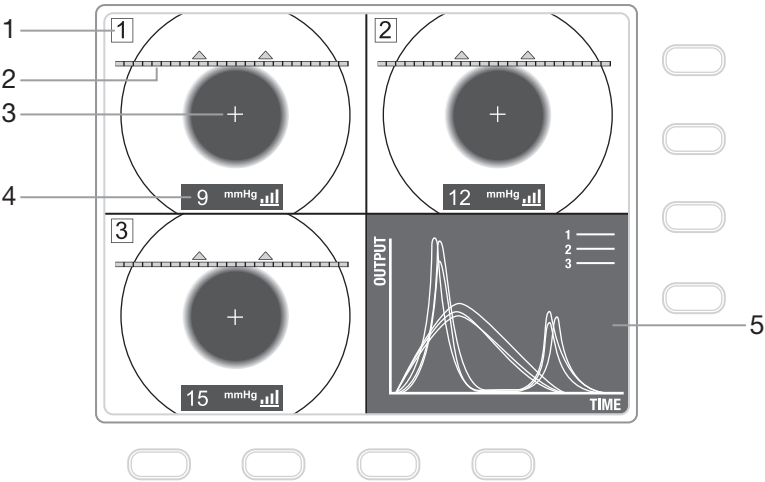


- | | | | |
|---|---|---|---|
| 1 | Measurement count display | 4 | Intraocular pressure measurement value
Intraocular pressure measurement value, reliability mark, received light level mark, and unit of measure are displayed. |
| 2 | Eyelid warning line
Perform the measurement while the eyelid and eyelashes are above the warning line. | 5 | Intraocular pressure measurement waveform
Waveform graph of pressure sensor signal and corneal distortion signal of measured values. |
| 3 | Cross mark
Displays the intraocular pressure measurement optical axis. | | |

Other Screens

USB Connection Check Screen

Screen for confirming that data is being input correctly from the USB device.

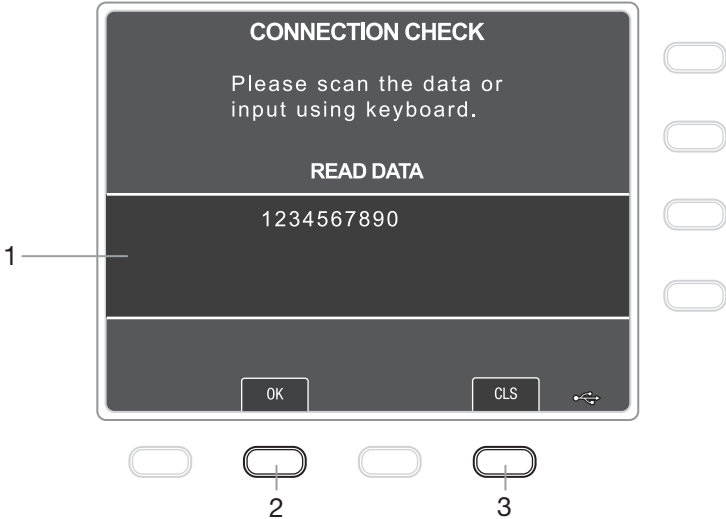

- 1


Data display field
The data input from the USB device is shown.
- 2

OK key
Press to close the USB connection check screen.
- 3

CLS key
Press to delete the data input from the USB device.

Print Message Input Screen

This screen is for inputting print messages.
Operations and the functions of the icons are the same as for the **PATIENT ID** screen.

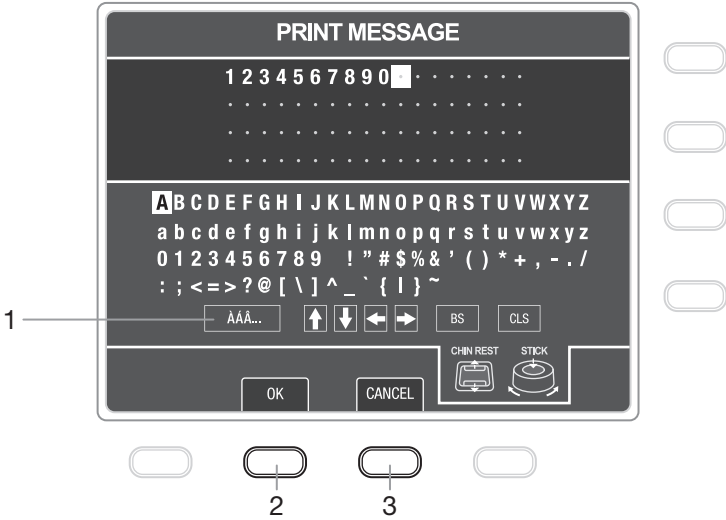

- 1


Input text switch command
AAA... : Switch the list of characters to Latin characters other than English.
ABC... : Switch the list of characters to English.
- 2

OK key
Press to save the print message and switch to the **MENU** screen.
- 3

CANCEL key
Press to not save the print message and switch to the **MENU** screen.

User Name Input Screen

Enter the user name associated with the computer that is connected via LAN in this screen.
Operations and the functions of the icons are the same as for the **PATIENT ID** screen.

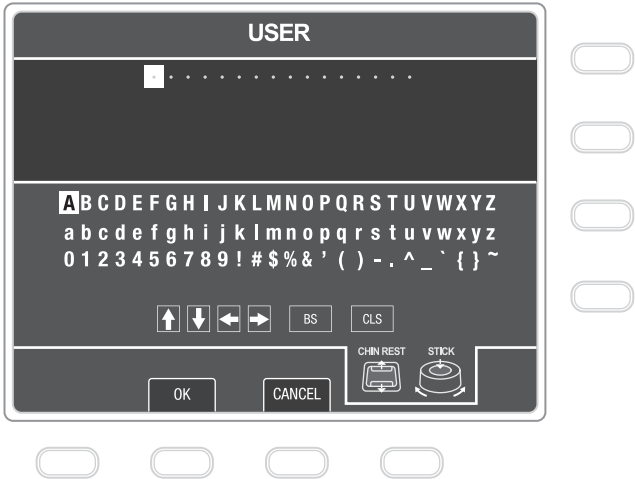


Password Input Screen

Enter the password for the computer that is connected via LAN in this screen.
Operations and the functions of the icons are the same as for the **PATIENT ID** screen.




Computer Name Input Screen

Enter the name associated with the computer that is connected via LAN in this screen. Input the name as text or an IPv4 format.
Operations and the functions of the icons are the same as for the **PATIENT ID** screen.

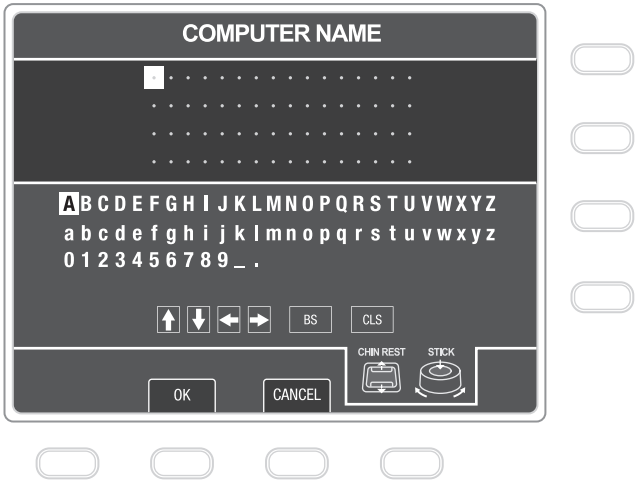

Output Folder Name Input Screen

Enter the name of a folder on the computer that is connected via LAN in this screen.
The slash is used as a separator so the setting can be deep in the hierarchy.
Operations and the functions of the icons are the same as for the **PATIENT ID** screen.

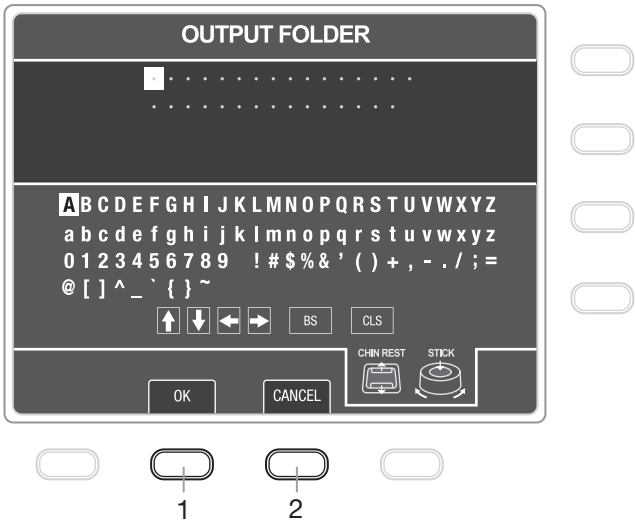


- 1

OK key
Press to save the output folder name and switch to the **MENU** screen.
- 2

CANCEL key
Press to not save the output folder name and switch to the **MENU** screen.

LCD Brightness Adjustment Screen


Adjust the LCD brightness in this screen.

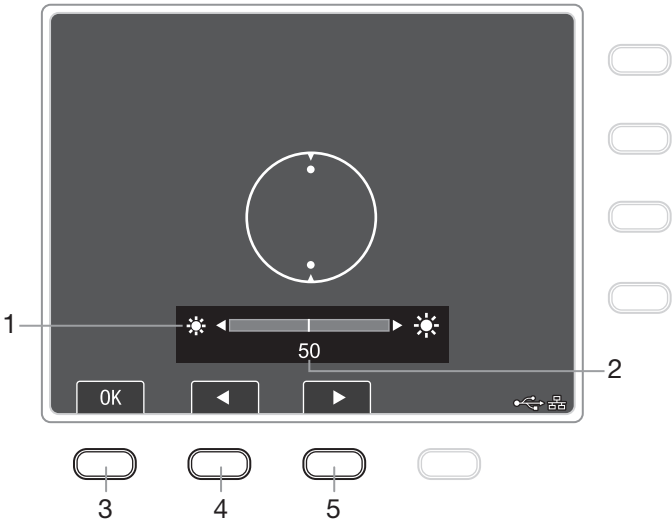

- | | | | |
|---|--|---|-------------------------------------|
| 1 | Brightness gauge | 3 | OK key |
| | Shows the current brightness of the LCD. The left side is 1 and the right side is 100. | | Switches to the MENU screen. |
| 2 | Brightness value | 4 | Brightness down key |
| | Shows the current brightness. Specify a value from 1 to 100. | | Reduce the brightness. |
| | | 5 | Brightness up key |
| | | | Increase the brightness. |

Description of Printouts

The following shows how measurement results are printed.

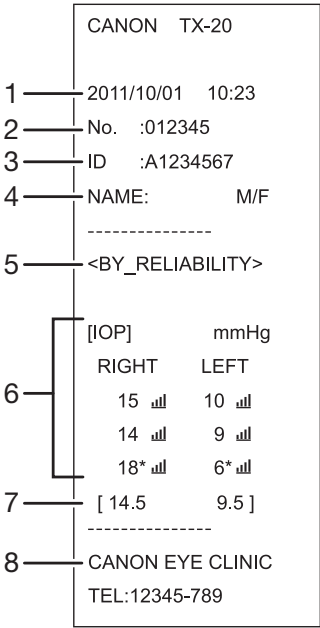

- | | | | |
|---|---|---|---|
| 1 | Date and time of measurement | 6 | Intraocular pressure measurement |
| 2 | Patient number | 7 | Intraocular pressure standard value |
| 3 | Patient ID | 8 | Message |
| 4 | Patient information entry field | | A message that has been input is printed. |
| 5 | Order of data | | |
| | <BY RELIABILITY>: In order of reliability | | |
| | <BY TIME>: In order of time | | |
| | <STANDARD ONLY>: Only the standard | | |

4 Basic Operation

This chapter explains the flow of operations and the basic measurement procedure.

Flow of Operations


The TX-20 has three alignment modes (full auto, auto, and manual).

Full auto: Alignment, measuring left and right eyes, and printing the results of the measurement are all done automatically.

Auto: Alignment and 1 measurement are done automatically. This mode is used to do measurements one eye at a time while checking the condition of the eye, such as for patients who cannot keep their eyes open continuously.

Manual: All operations are done manually. This mode is used when full auto or auto modes cannot be used, such as when the patient's pupil is eccentric or cornea is deformed.

The flow of operations is shown in the flow chart below. The items with gray backgrounds are performed automatically.

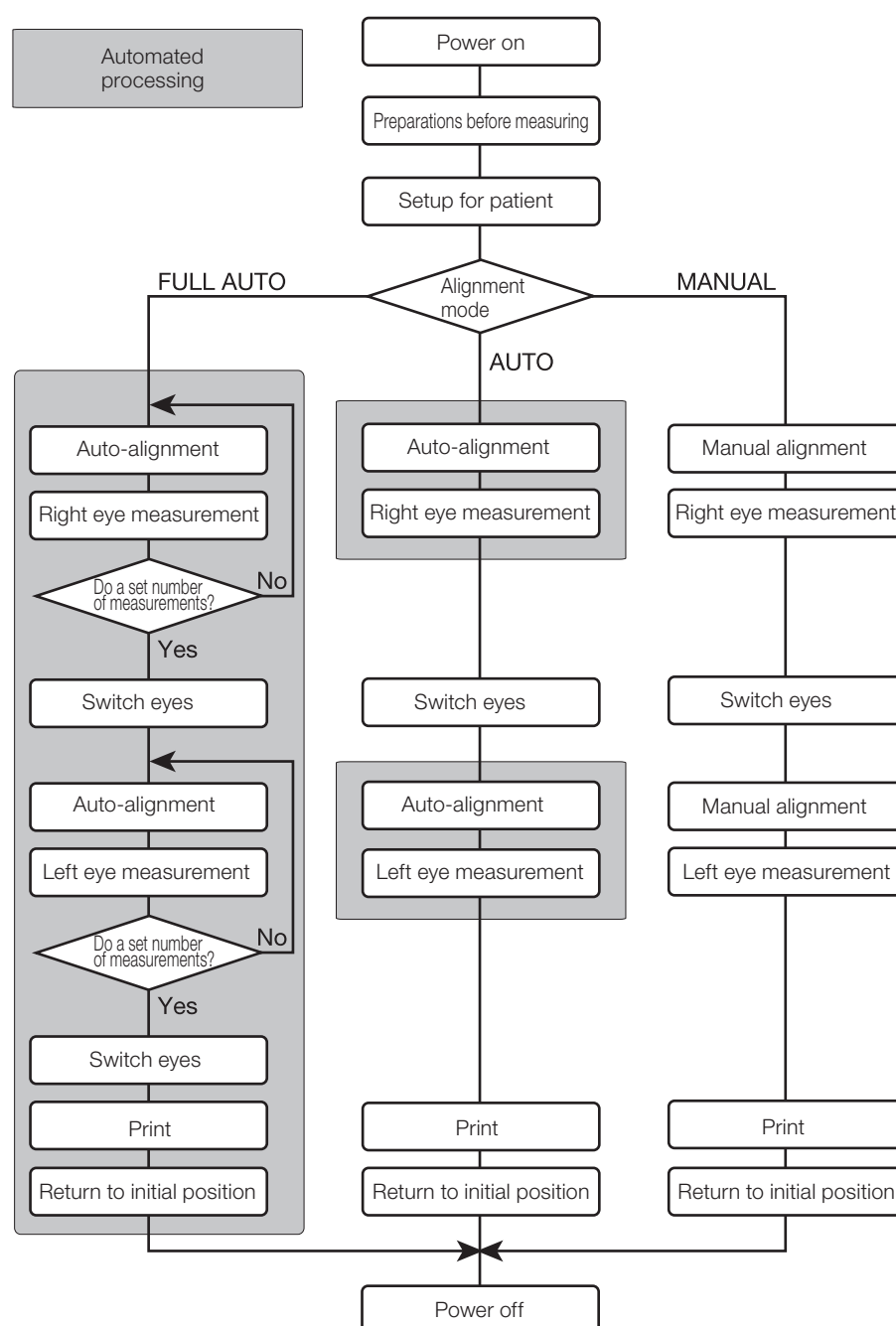

Connecting Cables


WARNING

- Take the following precautions to avoid risk of fire or electric shock:
 - Do not handle the power plug with wet hands.
 - Use only the power cord provided with the instrument.
- Connect the cables carefully so there are no problems operating the instrument.

CAUTION

- Always turn off the power and unplug the power cord before connecting the cables of any external devices to the TX-20.
- Use the computer and other equipment that conform to the system standard IEC 60601-1 or IEC 60950-1 for the Full Auto Tonometer TX-20. Be sure that the entire system conforms to IEC 60601-1. Be sure to also use an isolation transformer conforming to IEC 60601-1 when a computer conforming to IEC 60950-1 is used. Otherwise, electric shock may occur. Ask your sales representative or local Canon dealer to connect it to a computer.
- If the LAN connector is being used to connect to a network, connect an isolation transformer for networks between the TX-20 and network devices (HUB etc.) and between network devices. If this is not done, there is a risk of short circuits or malfunctioning devices. Ask your sales representative or local Canon dealer to connect it to a network.

1 Make sure that the power switch is off.

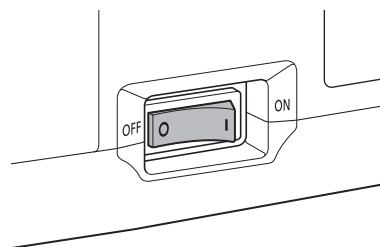


2 Connect a USB cable to input patient IDs with USB devices.



Remove the paper from the printer before turning the TX-20 on its side to attach cables. Otherwise the paper may jam.




- Always turn off the power and unplug the power cord before connecting the USB cable to the TX-20.
- For details on USB devices for inputting patient ID, contact your sales representative or local Canon dealer.

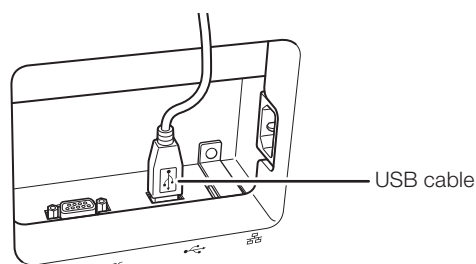


3 Connect an RS-232C cable or LAN cable to send measurement results to an external device.



- You cannot use the RS-232C and LAN connections at the same time. Use only one.
- Use a cross cable for the RS-232C cable.
- Ask your sales representative or local Canon dealer to set the connections for RS-232C and LAN.

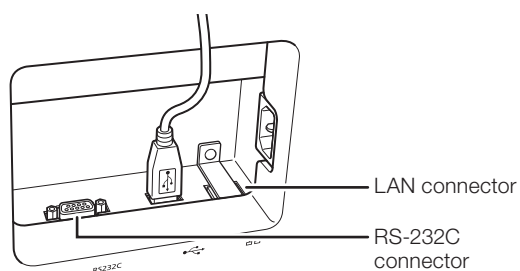

4 Connect the power cord.

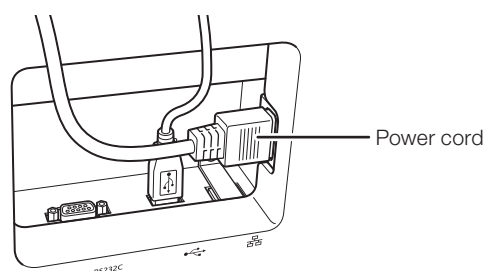

5 Insert the power plug all the way into an AC outlet.



Do not replace the paper unless the TX-20 is correctly installed in a flat location as specified. Changing the paper while the TX-20 is not level may cause a paper jam.

Preparations Prior to an Examination

CAUTION

Disinfect the forehead rest and replace the chin rest paper for each patient

To prevent the risk of infection, wipe the forehead rest with disinfectant ethanol for each patient. If the chin rest paper is not being used, be sure to disinfect the chin rest also. The forehead rest and chin rest may be corroded if a disinfectant other than those above is used. For details on how to disinfect, consult a specialist.

1 Remove the dust cover and the nozzle cap and check the printing paper.

Do not let dust get on the nozzle cap.


Add more paper for printing if there is not enough (see page 74).

2 Turn on the power.

CAUTION

The measurement unit moves during initialization-keep the patient's chin away from the chin rest and do not touch the measurement unit

There is a risk of injuring the patient.

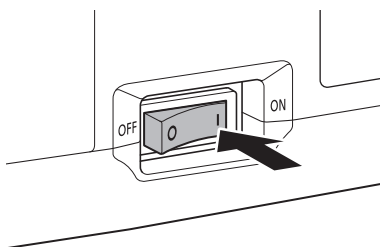

The instrument automatically performs an initialization sequence (internal check and positioning of measurement unit). After initialization, the measurement unit moves to the right eye side (initial position) and stops.




Initial position of measurement unit

You can set the initial position of the measurement unit to either the right or left eye (see page 61).

3 Adjust the tilt of the monitor.

Use both hands as shown in the illustration when adjusting the tilt of the monitor.

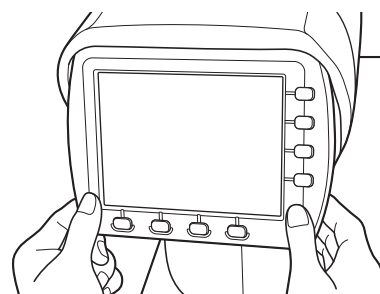


4 Disinfect the forehead rest and the chin rest.

If chin rest paper is being used, change the paper.

Patient Setup

- 1** Disinfect the forehead rest and the chin rest.
If chin rest paper is being used, change the paper.
- 2** Have the patient sit down.
- 3** Instruct patients to remove their glasses or contact lenses if they are wearing them.
 - Explain to the patient that intraocular pressure is measured by blowing air against the cornea.
 - Have the patient place their hand in front of the nozzle, and press the CHECK button, so they can feel the puff of air.
- 4** Input the patient ID.
Input the patient ID using a magnetic card reader, barcode reader, numeric keypad, or other USB device.
To input the patient ID on the TX-20 without using a USB device, press the **ID** key and input it from the **PATIENT ID** screen.
Go to step 5 if you are not going to input the patient ID.
- 5** Press the R/L button to decide which eye to measure.
The standard is to measure the right eye first. To measure the left eye, press the R/L button to switch the position of the measurement unit.
- 6** Place the patient's chin on the chin rest and forehead against the forehead rest.
Adjust the heights of the optical bench and chair so the patient is comfortable.
Be careful that the patient's chin does not protrude from the chin rest.
- 7** Move the chin rest with the CHIN REST button so that the patient's eye is aligned with the height adjustment mark.

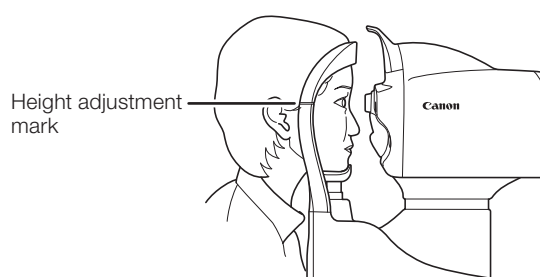


8 Press the LIMITER button to set the position where the nozzle stops.

CAUTION


Be sure to set the position where the nozzle stops for each patient

Look from the side of the patient when setting the stop position to see that the patient's forehead and chin are firmly against the forehead rest and chin rest. If the stop position is not set correctly, the nozzle may contact the patient's eye and result in injury.

- 1) Check that the LIMITER OFF mark is displayed.

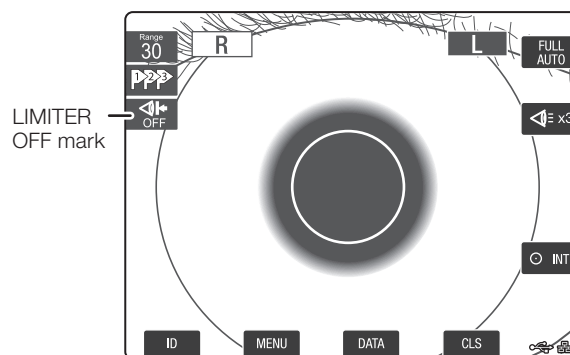


- 2) Slowly move the operation lever so the nozzle approaches the eye being examined.
Look at the patient from the side to make sure that the patient's forehead and chin are in secure contact with the forehead rest and chin rest.
- 3) Set the distance between the nozzle and the patient's eye to approximately 8 to 10 mm, and press the LIMITER button.
The LIMITER OFF mark disappears and the nozzle stop position mark is displayed.

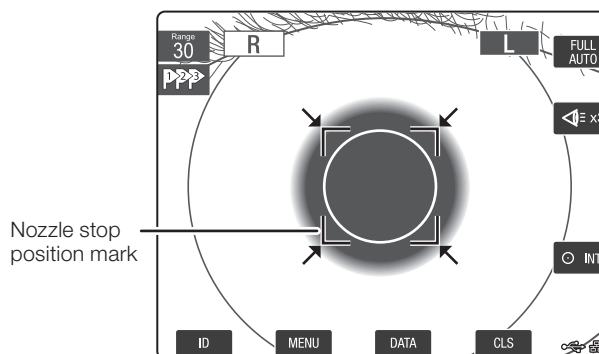

Stop doing measurements if the LIMITER OFF mark does not disappear even if the LIMITER button is pressed

Stop using the instrument, and contact your sales representative or local Canon dealer.

- 4) Move the operation lever to check the following:
- Pull the operation lever slightly backward so the nozzle stop position mark disappears.
 - Tilt the operation lever forward and the nozzle stop position mark appears. The nozzle should not approach the eye being examined closer than the stop position.

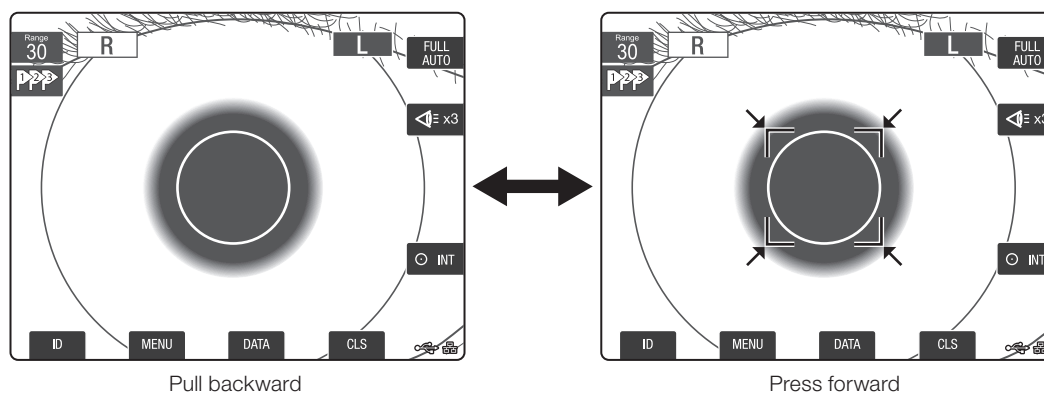


9 Have the patient stare at the eye fixation lamp.

Have the patient stare at the green internal eye fixation lamp inside the nozzle.



Press the eye fixation lamp key to switch to the external eye fixation lamp if the patient cannot see the internal eye fixation lamp, if it is difficult to see, or if the patient's line of sight is unsteady.

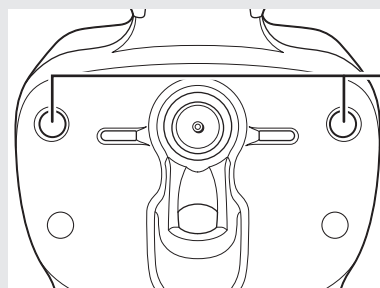

External eye
fixation lamp

Measuring

This section explains the full auto-alignment mode, which is most often used for normal measuring. The number of measurements is assumed to be set at three.




To get a highly reliable measurement

We recommend that the number of measurements be set to three in full auto-alignment mode. We also recommend to perform more than three measurements in the auto- or manual-alignment mode to get a stable measured values.

4 Basic Operation

- 1 Use the operation lever and the measurement unit vertical movement ring to display the patient's pupil in the monitor.

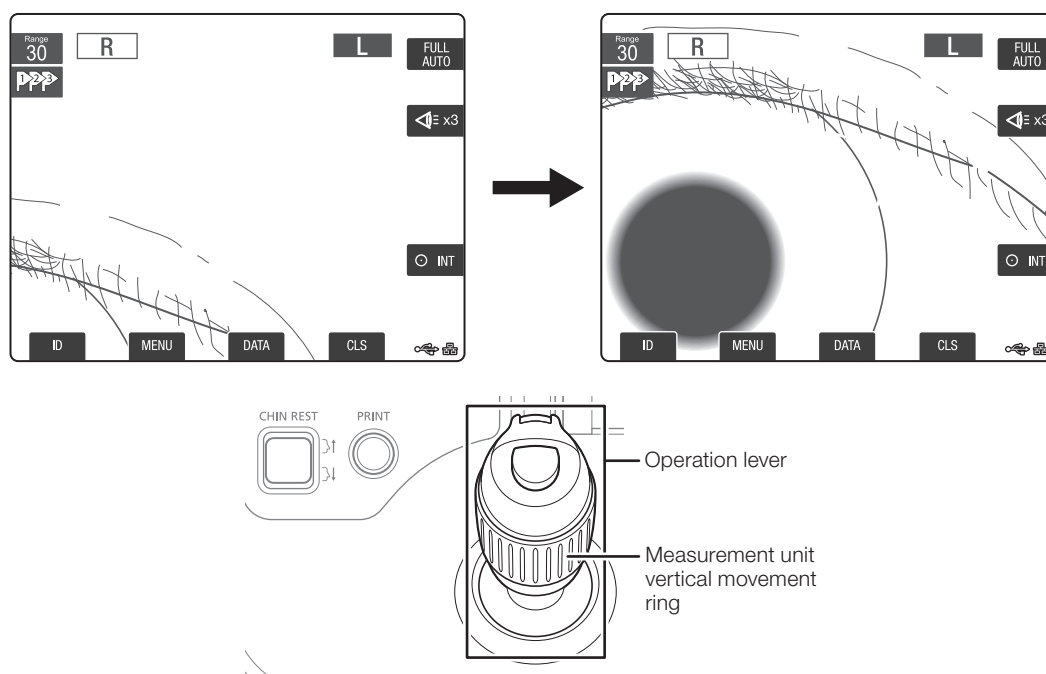

During full auto/auto-alignment

Show all of the pupil in the monitor.

During manual alignment

Have the alignment circle overlap close to the center of the pupil. Also, align the corneal light point on the same straight line as the vertical alignment guide (see page 45).



How to use the operation lever

- Left tilt: The measurement unit moves to the left a little according to the left tilt. The motion is greater and speed increases as the tilt surpasses the maximum tilt position (the buzzer sounds briefly).
- Right tilt: The measurement unit moves to the right a little according to the right tilt. The motion is greater and speed increases as the tilt surpasses the maximum tilt position (the buzzer sounds briefly).
- Forward tilt: The measurement unit moves toward the patient a little according to the forward tilt. The motion is greater as the tilt surpasses the maximum tilt position (the buzzer sounds briefly).
- Backward tilt: The measurement unit moves toward the examiner a little according to the backward tilt. The motion is greater as the tilt surpasses the maximum tilt position (the buzzer sounds briefly).
- Rotate right: Measurement unit moves up.
- Rotate left: Measurement unit moves down.



To adjust the angle of the operation lever without moving the measurement unit

Press the alignment off button while adjusting the angle of the operation lever. While the button is pressed, the measurement unit does not move even if the operation lever is tilted.


2 Have the patient stare at the eye fixation lamp.

Ask patients if the green internal eye fixation lamp can be seen and instruct them as follows:

- Have them blink several times before measurement.
- Have them open their eyes very wide during measurement.
Help them to keep their eyelid open if the movement is difficult for them. Be careful not to apply pressure to the eyeball while holding the eyelid open. It is not possible to do measurements if there is pressure on the eyeball.
- Have them not blink during measurement.
- Have them stare at the eye fixation lamp.
Switch to the external eye fixation lamp if it is difficult to see the internal eye fixation lamp or if the patient's line of sight is unsteady (see page 43).

3 Press the measurement button.

- 1) The measurement unit moves automatically to align with the position of the patient's eye.

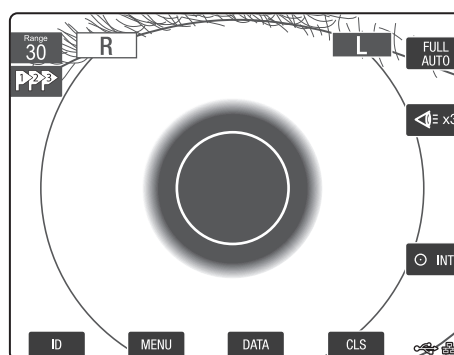

To stop measurement

- Press any button on the operation panel.
- Press any soft key on the monitor panel.
- Press the measurement button or the alignment off button.
- Tilt the operation lever.
- Rotate the measurement unit vertical movement ring.



To do a measurement before alignment is finished

- Press the measurement button.

Press the measurement button when the corneal light point can be seen clearly. You can do measurements even during auto-alignment.

However, the positioning alignment may not be accurate. Note that the measurement may not be accurate when measuring this way.



To restart measurement

- Press the measurement button.

If you are doing a set number of measurements, when the measurement button is pressed the eye to which the measurement unit is positioned is measured once.



If a message appears

Check "7 Troubleshooting" (see page 77) and perform the remedy indicated in the message.

- 2) Three consecutive measurements are performed.
When the position is aligned, three measurements are done consecutively.
- 3) The measurement unit moves to the other eye.
The set number of measurements is done and the measurement unit moves to the other eye automatically.



Except for the full auto-alignment mode, switching between the left and right eye is done manually.

- 4) Alignment is done with the left eye; three measurements are done consecutively.
- 5) The measurement results are printed (see page 36).



Sending data to an external device

If you have set the instrument to send measurements to an external device, the measurement results are output at this point.

If all the measurements are erroneous

The measurement values are neither printed out nor output via RS-232C. However, the erroneous data is output via LAN.

Stopping printing or data transfer in progress

To stop printing or data transfer in progress, press a button on the operation panel or a soft key on the monitor panel, measurement button, or alignment off button.

- If you stop printing, the paper is not cut automatically.
- If you stop transferring data, the file after the file being transferred is not sent.



The thermally sensitive paper used for printing may degrade over time. Make a copy of measurement results to store them for a long period.

- 6) The measurement unit moves to the right eye.
The **Measurement** screen changes to the default state. The measured value is cleared and the LIMITER OFF mark is displayed.
However, you can check the measured value on the **DATA** screen. When the measurement button is pressed for the next measurement, the measured value on the **DATA** screen is also cleared.
- 7) Repeat the process from "Setup for Patient" to continue measurements with a different patient.

Finishing an Examination

- 1 Turn off the power switch.
- 2 Clean the nozzle.
- 3 Disinfect the forehead rest and the chin rest.
If chin rest paper is being used, change the paper.
- 4 Put the cap on the nozzle and cover the equipment with the dust cover.



To pack and transport

Turn on the power in packaging mode (see page 76), let the measurement unit move to the packaging and transport position, and then turn off the power.

Remove the power cable and remove the cable from any external devices.

Power Saving Mode

If no operations are performed for a set period, the TX-20:

- Turns off the backlight of the LCD monitor.
- Turns off all LED lamps, except the power lamp.
- Flashes the power lamp.

Recovering from the Power Saving Mode

Press a button on the operation panel, a soft key on the monitor panel, or the measurement button, or tilt the operation lever to recover from the power saving mode.

Either the **Measurement** screen or the **PATIENT ID** screen is displayed at recovery, depending on the [Redisplay In Scn] setting in the Patient ID menu (see page 68).



The status of the TX-20 is the following when recovering from power saving mode:

Do the following operations if necessary.


- Nozzle stop position: The setting value is released. Follow the procedure and set the position again.
- Eye fixation lamp: The internal eye fixation lamp is selected. If the external eye fixation lamp was selected, reselect it.
- Patient No.: Counts up the patient numbers. (It depends on the settings.) Set the number again if required.
- Measurement screen: The displayed value is cleared. Check the patient ID or measured value on the DATA screen or perform the measurement again.

5 Types of Measurements

This chapter explains how to perform the settings on the **Measurement** screen, **PATIENT ID** screen, and **MENU** screen.

Measurement Screen

You can set the alignment mode, number of measurements, and the eye fixation lamp in the measurement mode.

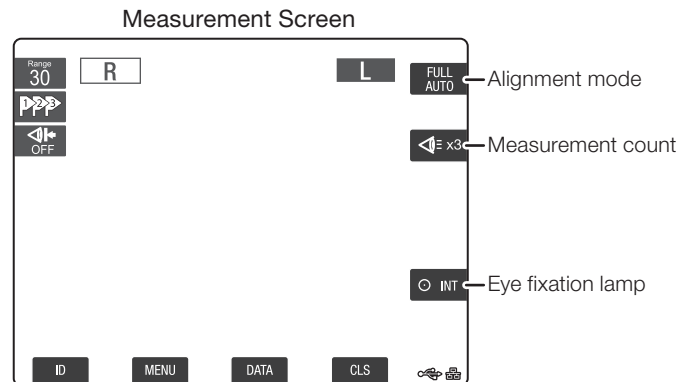


Changing the Alignment Mode

The TX-20 lets you choose between the full auto, auto, and manual alignment modes.

- **Full auto-alignment mode (FULL AUTO)**
Press the measurement button to do alignment, measurement (1 to 3 times), switching eyes, alignment, measurement (1 to 3 times), and output of measurement results automatically.
In the **MENU** screen you can set whether or not to automatically output the measurement results.
- **Auto-alignment mode (AUTO)**
Press the measurement button to do alignment and one measurement automatically.
- **Manual-alignment mode (MANUAL)**
All operations are done manually.

- 1 Press the alignment mode key.
Every time the key is pressed, the alignment mode changes.

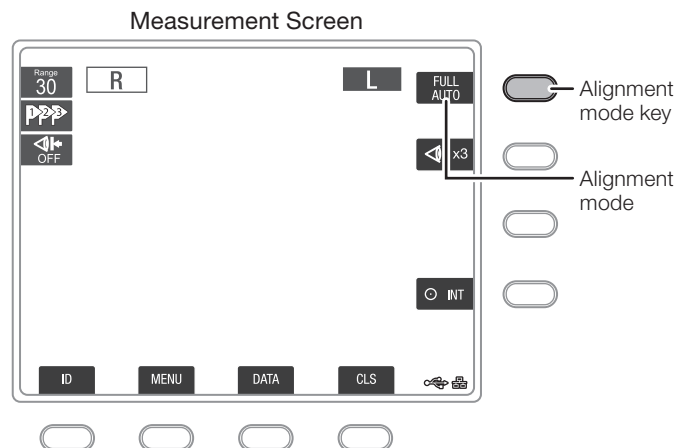


Changing the Number of Measurements

You can set the number of measurements done in full auto-alignment mode.

- One measurement ($\triangleleft=x1$): One automatic measurement is done.
- Two measurements ($\triangleleft=x2$): Two automatic measurements are done.
- Three measurements ($\triangleleft=x3$): Three automatic measurements are done.

- 1 Press the measurement count key.
Every time the key is pressed, the number of measurements changes.

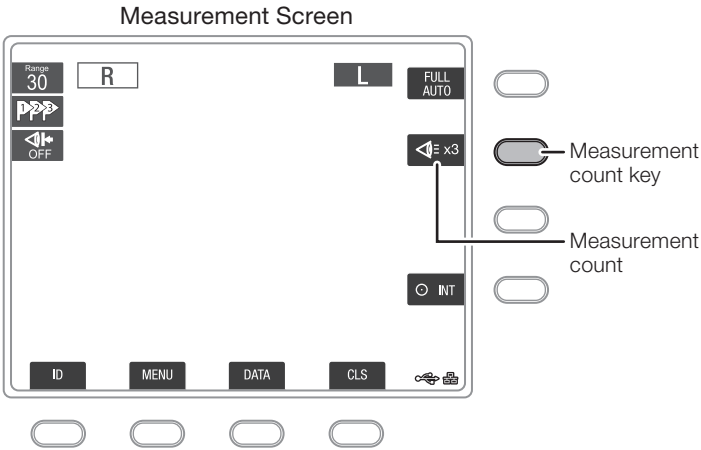


To get a highly reliable measurement


We recommend that the number of measurements be set to three in full auto-alignment mode.

We also recommend performing more than three measurements in the auto- or manual-alignment mode to get a stable measured values.


Changing the Eye Fixation Lamp

You can select either the internal or external eye fixation lamp.

Internal eye fixation lamp (): Located inside the nozzle. Have the patient stare to fix his/her line of sight. Normally, the internal eye fixation lamp is used.


External eye fixation lamp (): Located outside the nozzle. The external eye fixation lamp stabilizes the line of site for patients who cannot see the internal eye fixation lamp or whose line of site is not stabilized by the internal eye fixation lamp.

- 1
- Press the eye fixation lamp key.
- Every time the key is pressed, the eye fixation lamp changes.

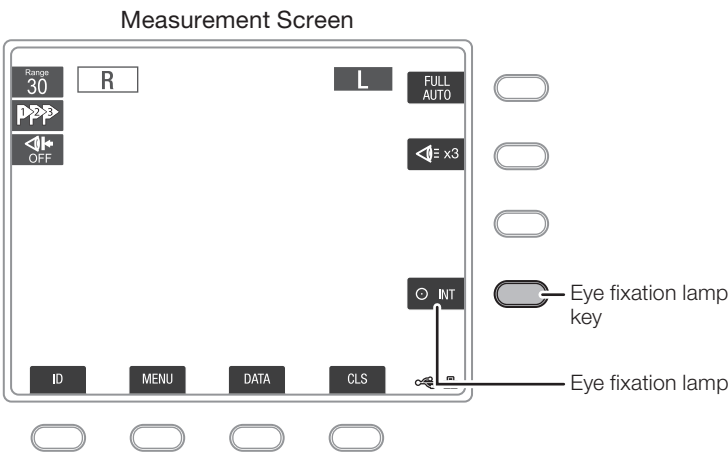


PATIENT ID Screen

Set the patient ID in the **PATIENT ID** screen. The patient ID is a code to identify patients; you can set up to 64 alphanumeric characters and symbols. The patient ID can be displayed to appear on the **Measurement** screen, **PATIENT ID** screen, printed data, and data output to external devices. You can input patient IDs from a bar code reader or magnetic card reader that is connected via the USB interface (Human Interface Device (HID) class).



Patient ID on the Measurement screen


Six digits of the patient ID can be shown on the **Measurement** screen. If a patient ID is longer than 6 digits, "..." appears at the end of the ID. Confirm the entire patient ID in the **PATIENT ID** screen.

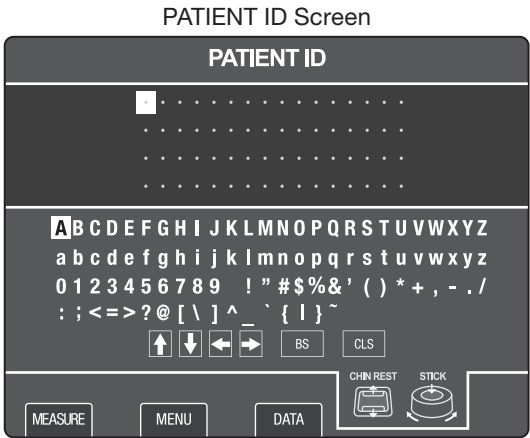

Input Patient ID in the PATIENT ID Screen

This section explains how to input the patient ID in the **PATIENT ID** screen.

PATIENT ID Screen Display and Functions of Buttons

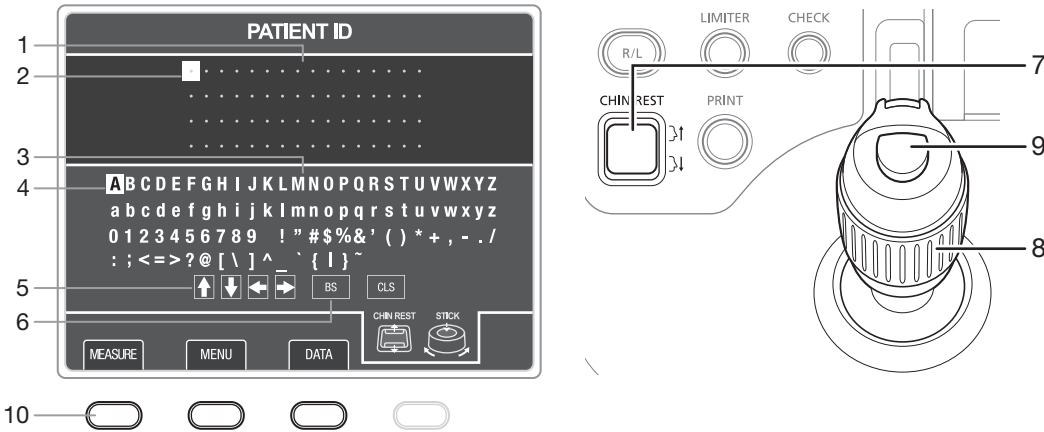


- 1 Patient ID area
Displays the patient ID that was input.
- 2 ID cursor
Displays the position immediately after a character is input. Also, selects the targeted character when deleting a patient ID.
- 3 List of characters
List of characters that can be input as a patient ID. Specify from upper and lower case letters, numbers, and symbols.
- 4 Selection cursor
Selects characters or commands to use.
- 5 ID cursor movement commands
Moves the ID cursor up, down, left, or right one space.
- 6 ID delete command
Deletes the patient ID.
BS : Deletes the character that is before the character specified by the ID cursor.
However, if the ID cursor is specifying the first character, it deletes that character.
CLS : Deletes all the characters in the patient ID.
- 7 CHIN REST button
Moves the selection cursor up and down.
- 8 Measurement unit vertical movement ring
Moves the selection cursor left and right.
Clockwise: Right
Counterclockwise: Left
- 9 Measurement button
Inputs the character specified by the selection cursor for a patient ID.
- 10 **MEASURE** key/**MENU** key/**DATA** key
Moves to other screens.
MEASURE : **Measurement** screen
MENU : **MENU** screen
DATA : **DATA** screen

Basic Operations on the PATIENT ID Screen


This section explains the basic procedures for the patient ID settings.

- Steps 2 to 4: Input text.
- Steps 5 to 6: Delete text.
- Steps 7 to 8: Add text to existing text.
- Step 9: Delete the entire ID.

- 1 Display the **PATIENT ID** screen.
Press the **ID** key.
The ID cursor is at the end of the patient ID on the right.
The selection cursor is at the top of the list of characters.

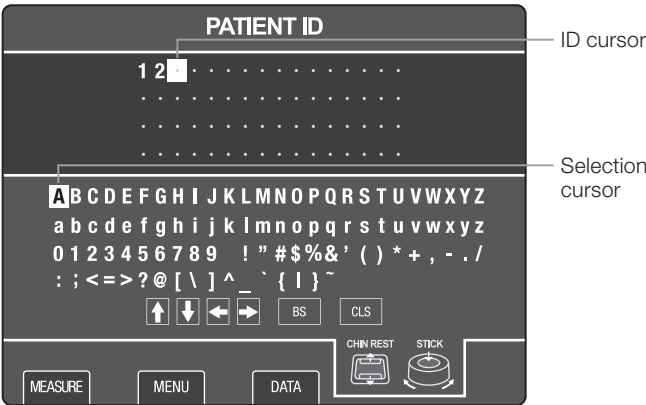


- 2 Decide on the position to input text.
Move the ID cursor to the right of where you want to input a character.
The procedure in steps 2 to 4 is explained using an example of inputting a "C" at the end of patient ID "12".
In this example, the cursor is already to the right of the "2", so there is no operation.

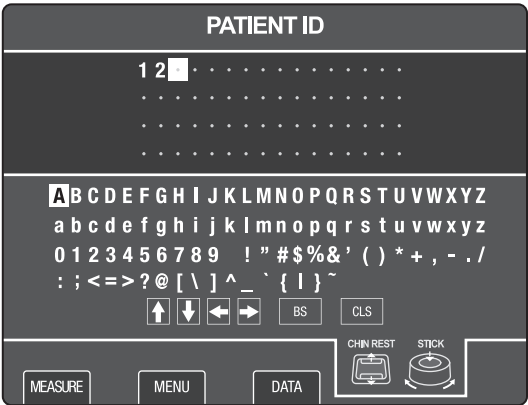


To add text to the front
Move the ID cursor to the first character.


- 3
- Select a character to input.
- Use the CHIN REST button and the measurement unit vertical movement ring to move the selection cursor to the position of the [C] in the character list.

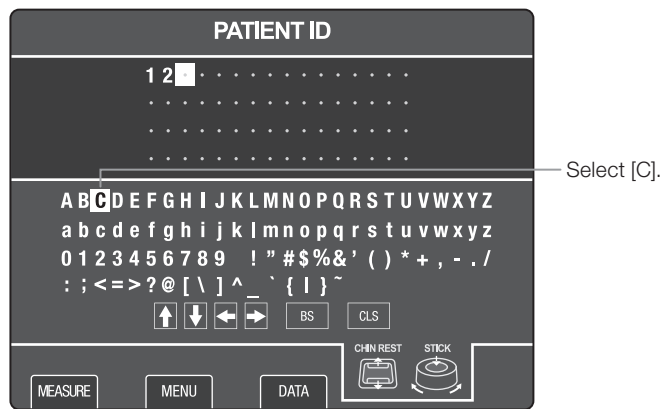

- 4
- Input the selected character.
- While [C] is selected with the selection cursor, press the measurement button one time.
- The "C" appears to the right of the "2".

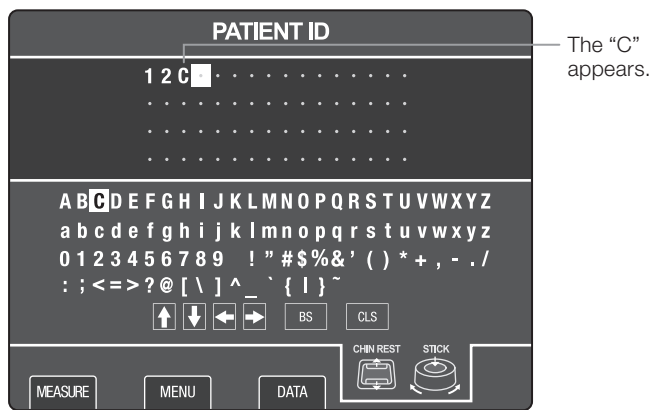


- 5
- Select the place to the right of the character to delete.
- Move the ID cursor to the right of the character to delete.
- Use the CHIN REST button and the measurement unit vertical movement ring to move the selection cursor to the position of [←] (ID cursor movement command). Press the measurement button one time to move the cursor one space to the left.
- The procedure in steps 5 to 6 is explained using an example of deleting a "2" from the patient ID "12C". In this example, the cursor moves to the "C" to the right of the "2".

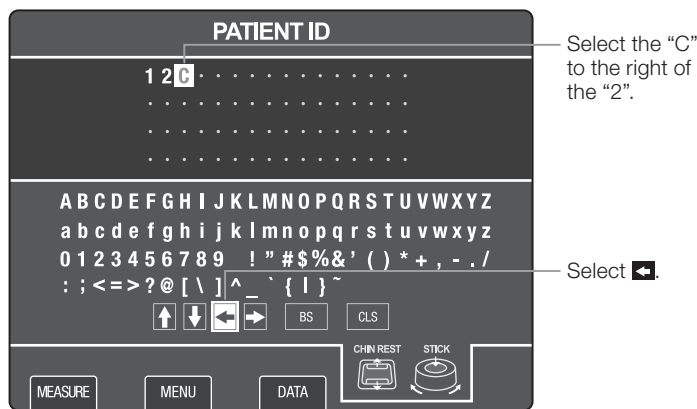


6 Delete a character.


Use the CHIN REST button and the measurement unit vertical movement ring to move the selection cursor to the position of **BS** (ID delete command).
Press the measurement button one time to delete “2” and move “C” one space to the left.

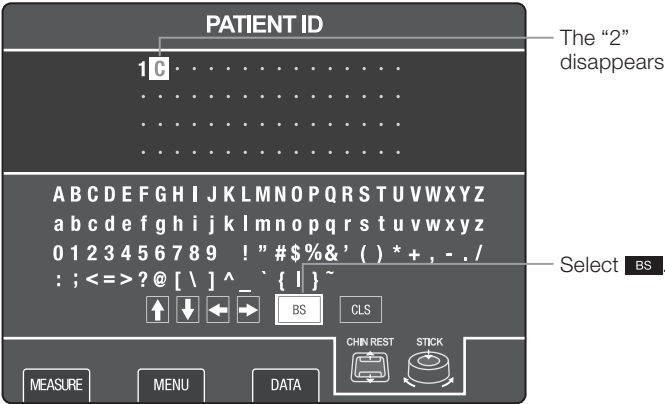

7 Decide on the position to add text.


Move the ID cursor to the character to the right of the place to add text.
The procedure in steps 7 to 8 is explained using an example of adding an “@” between the “1” and “C”.
In this example, the cursor is already at the position of the “C”, so there is no operation.

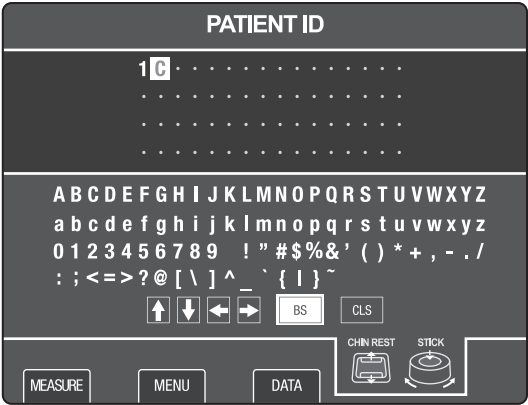


To add text to the front
Move the ID cursor to the first character.


8 Add a character.

Use the CHIN REST button and the measurement unit vertical movement ring to select **[@]**. Press the measurement button one time to add the “@” in front of the “C”.

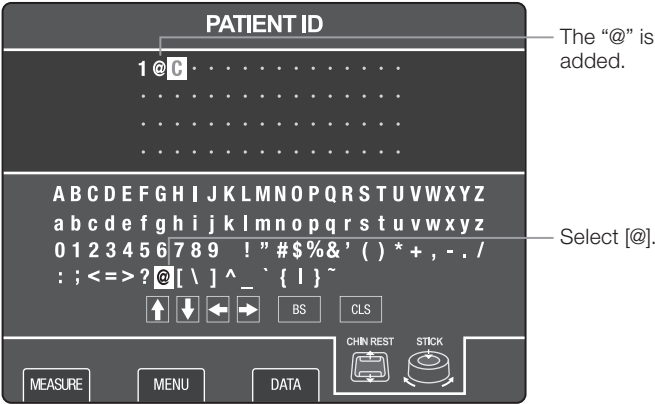

9 Delete the entire patient ID.

Use the CHIN REST button and the measurement unit vertical movement ring to move the selection cursor to the position of **CLS** (ID delete command).
Pressing the measurement button opens a message screen to confirm that you want to delete the entire patient ID.
Select [OK] to delete. All of the “1@C” patient ID is deleted.
Select [CANCEL] to not delete. The “1@C” patient ID is kept as is.

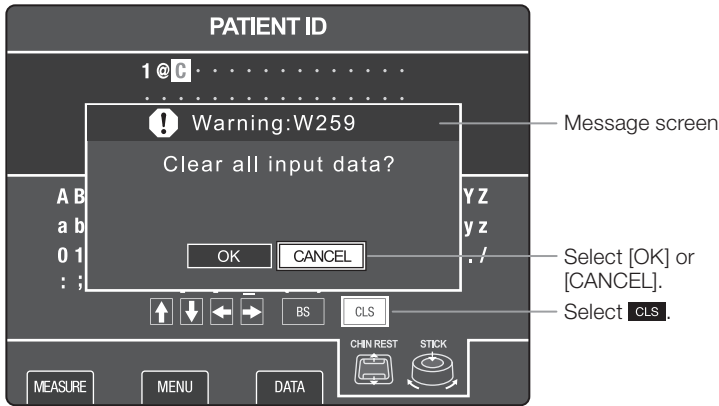

Input Patient ID from an External Device

Before starting measurements, input the patient ID according to the operating procedure for the input device.


For information about the TX-20 and connections with a USB device, settings and usage, contact your sales representative or local Canon dealer.

- The keys on a USB device that can be input are shown below.
- Numbers
 - Letters (upper and lower case)
 - Symbols (excluding “#”, “%”, “_”, and backslash)
 - Space, BackSpace, Delete, and arrow keys
 - NumLock and Enter

MENU Screen

In the **MENU** screen, you can set the measurements, printing, data input and output, patient ID, patient No., date & time, messages, and the brightness of the monitor.

MENU Screen Display and Functions of Buttons

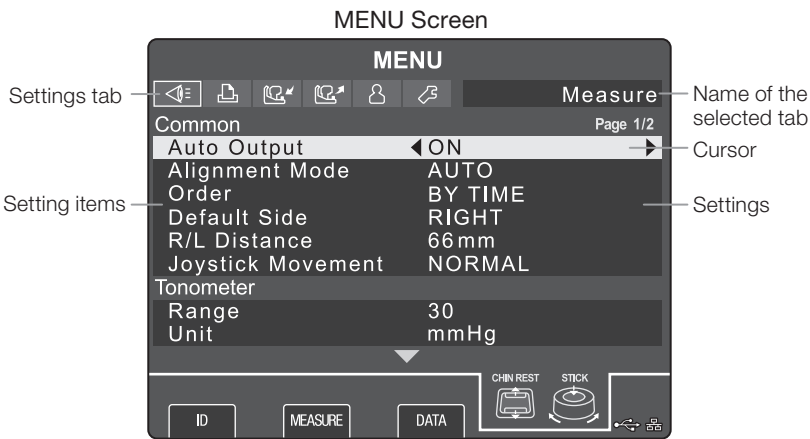

- **Settings tab**
The items to set are grouped according to their purpose. Selected items are light in orange.
 - Measure tab: Settings for intraocular pressure measurements
 - Print tab: Settings for printing with the built-in printer
 - Input tab: Settings for inputting data from an external device
 - Output tab: Settings for outputting data to an external device
 - Patient tab: Settings for patient ID and patient No.
 - Other tab: Settings for time, power saving mode, the LCD brightness, and shows information about the firmware.
- **Setting items**
The setting items for each setting tab are displayed.
- **Settings**
The settings for each setting item are displayed.
- **Cursor**
Shows the position of the selected setting value, setting item, and tab. They light orange when selected.
- **CHIN REST button**
Moves the cursor up and down.
- **Measurement unit vertical movement ring**
 - Moves the cursor left and right.
 - Clockwise: Right
 - Counterclockwise: Left
 - Changes the settings and their values.

5 Types of Measurements


- **Measurement button**
 1. Moves the cursor between digits in the patient No., IP address, subnet mask, default gateway, DNS server, time and date settings. The cursor moves to the right every time it is pressed.
 2. Execute the setting.
 - [Input Message]: Move to the print message input screen.
 - [Adjust Brightness]: Move to the LCD brightness adjustment screen.
 - [Reset]: Set the patient No. to "1".
 - [Connection Check]: Check the connection with external devices.
- **Alignment off button**


Moves the cursor between digits in the patient No., IP address, subnet mask, default gateway, and DNS server settings. The cursor moves to the left every time it is pressed.


Basic Operations on MENU Screen

This section explains how to set a patient No. You can set other items with the same operations. This procedure shows how to change the patient No. from "000000" to "000012".

- 1 Display the **MENU** screen.


Press the **MENU** key in the screens.
- 2 Select the  Patient tab.

Use the CHIN REST button and the measurement unit vertical movement ring to select  Patient tab.

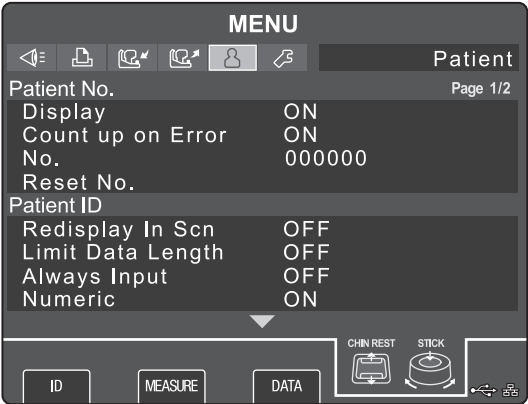


- 3 Select an item.

Use the CHIN REST button to select [Patient No.] > [No.].

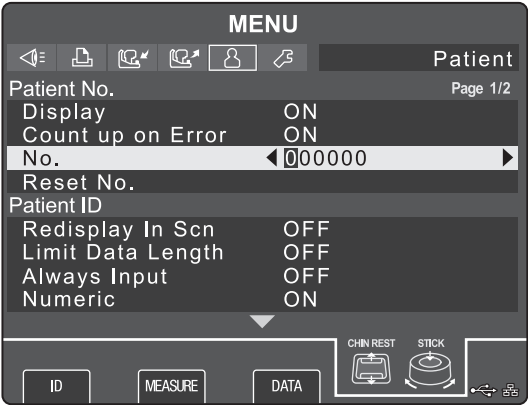


4 Select a digit.

Press the measurement button four times. Move the cursor four digits to the right and select a digit to set a number.
If you go too far, press the alignment off button to move the cursor to the left.

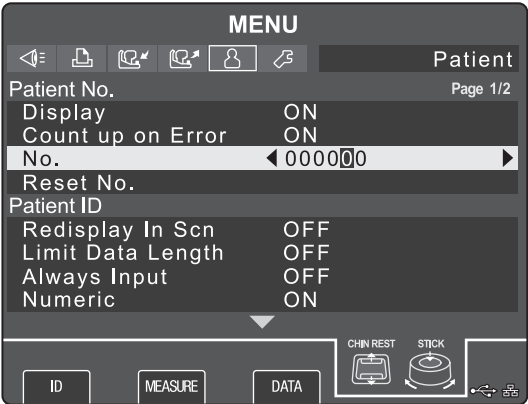


5 Input a value.

Rotate the measurement unit vertical movement ring left and right to set [1] as the selected digit.

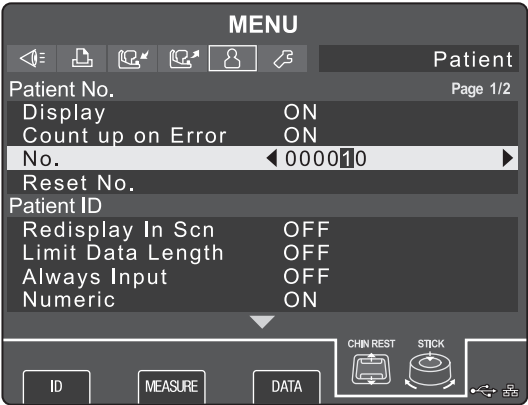


6 Select a digit.

Press the measurement button one time. Move the cursor one digit to the right and select a digit to set a number.
If you go too far, press the alignment off button to move the cursor to the left.

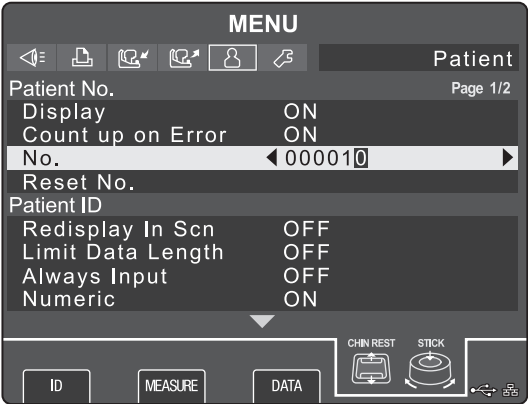

5 Types of Measurements

- 7
- Input a value.
Rotate the measurement unit vertical movement ring left and right to set [2] as the selected digit.

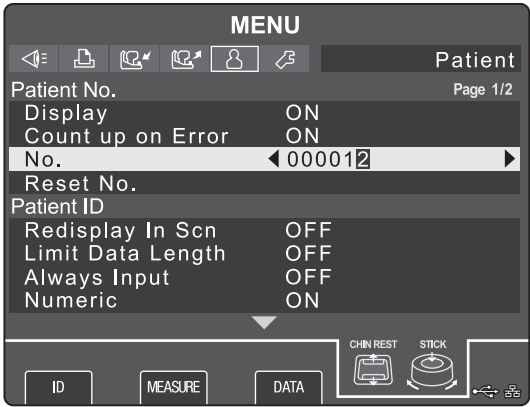

- 8
- Close the **MENU** screen.
Press the **ID** key, **MEASURE** key, or **DATA** key to open a different screen.

Common Settings in Measurement Modes

In the common settings menu, you can set auto-data output, alignment mode initial settings, measurement results display order, measurement unit recovery position, distance moved when switching eyes, and operation settings of the operation lever.

How to Enter the Common Settings Menu

MENU screen >  Measure tab > [Common]

Setting items	Description	Factory settings
Auto Output	<p>Select whether to print (or output data) results automatically after measuring both eyes, in full auto-alignment mode.</p> <p>[ON]: Results are printed automatically. [OFF]: Results are not printed automatically.</p> <p>If the [Output] setting is [OFF], printing is not done even if this setting is [ON] (see page 63). If the [Interface] setting is [OFF], data is not output even if this setting is [ON] (see page 64).</p>	ON
Alignment Mode	<p>Select an alignment mode after initialization.</p> <p>The main unit performs initialization in the following situations.</p> <ul style="list-style-type: none"> • The power is turned on. • The PRINT button is pressed. • The CLS key is pressed. • The TX-20 recovers from power saving mode. <p>[FULL AUTO]: Full auto mode [AUTO]: Auto mode</p>	FULL AUTO
Order	<p>Select the order in which measurement results are shown.</p> <p>[BY TIME]: Show in order of time. [BY RELIABILITY]: Show in order of reliability.</p>	BY TIME
Default Side	<p>Set the right eye side or left eye side as the recovery position for the measurement unit.</p> <p>The measurement unit moves to the recovery position of the side set in the following situations.</p> <ul style="list-style-type: none"> • Power on • After pressing the PRINT button • After pressing the CLS key <p>[RIGHT]: Right eye [LEFT]: Left eye</p>	RIGHT
R/L Distance	<p>Set the distance moved when switching between left and right eyes.</p> <p>[58 to 66]: Distance moved (mm)</p>	62
Joystick Movement	<p>Set whether to put the operation lever into proximity limiter mode.</p> <p>In the proximity limiter mode, coarse adjustment motions are prevented when the measurement unit moves within 2 mm of the alignment end position. Beyond that, fine adjustment motion takes place until the unit reaches the position.</p> <p>[NORMAL]: Normal mode [SAFETY]: Proximity limiter mode</p>	NORMAL

Settings for Intraocular Pressure Measurements

In the intraocular pressure settings menu, you can configure the settings of the intraocular pressure measurement range switching, units of intraocular pressure measurements, display of received light level, warning display threshold for intraocular pressure, display of error screen, and display of warning screen/waveform.

How to Enter the Intraocular Pressure Measurement Menu

MENU screen >  Measure tab > [Tonometer]

Setting items	Description	Factory settings
Range	Switch the range of intraocular pressure that can be measured. [30]: Measures a range from 0 to 30 mmHg. If the pressure is over 30 mmHg, it switches to 60 mmHg. [60]: Measures a range from 25 to 60 mmHg. If the pressure is under 25 mmHg, it switches to 30 mmHg.	30
Unit	Set the units for intraocular pressure measurements. [mmHg]: Units of measurement are shown in mmHg. [hPa]: Units of measurement are shown in hPa.	mmHg
Level Meter	Set to enable showing the received light level. [ON]: Show the received light level. [OFF]: Do not show the received light level.	ON
Upper IOP Warning	Set the upper pressure limit at which measurement is stopped (only enabled in full auto mode). [17 to 25 mmHg]: Upper limit of intraocular pressure at which to stop measurement [23 to 33 hPa]: Upper limit of intraocular pressure at which to stop measurement [OFF]: Do not show message or stop measurement.	18
Lower IOP Warning	Set the lower pressure limit at which measurement is stopped (only enabled in full auto mode). [1 to 9 mmHg]: Lower limit of intraocular pressure at which to stop measurement [1 to 12 hPa]: Lower limit of intraocular pressure at which to stop measurement [OFF]: Do not show message or stop measurement.	OFF
Error Display	Set whether to show an anterior segment image when a measurement error occurs. [ON]: Shows image. [OFF]: Does not show image.	ON
Warning Display	Set whether to show three anterior segment images and waveform when measurements are inconsistent. [ON]: Shows images and waveform. [OFF]: Does not show image and waveform.	ON

Settings for Printing

You can enable or disable printing, and set the text size, the density of print, and the type of print data, as well as whether to print a title, Patient ID, Patient No., and a message, in the print menu.

How to Enter the Print Menu

MENU screen >  Print tab > [Print]

Setting items	Description	Factory settings
Output	Set whether to print the measurement results with the built-in printer. [ON]: Prints. [OFF]: Does not print.	ON
Reduced	Set whether to reduce distance between lines of print. [ON]: Reduce the distance between lines. [OFF]: Do not reduce the distance between lines.	OFF
Density	Set the print density. [1 to 3]: Print density (higher numbers are more dense)	2
Data	Set which measurements to print. [STANDARD]: Only the standard [ALL]: Measurements for each of the set measurements and the standard Up to 10 sets of measurements are printed. Pressing and holding the PRINT button for two seconds or more prints [ALL], regardless of the setting.	ALL
Title	Set whether to print the names of the manufacturer and machine. [ON]: Prints. [OFF]: Does not print.	ON
Patient No.	Set whether to print the Patient No. [ON]: Prints. [OFF]: Does not print.	ON
Patient ID	Set whether to print the Patient ID. [ON]: Prints. [OFF]: Does not print.	ON
Message	Set whether to print a message. [ON]: Prints. [OFF]: Does not print.	OFF
Ext. Character	Set the extension character used in the message. [LATIN]: Latin Alphabet No. 1 / ISO 8859-1 standard [JIS]: Katakana	LATIN
Input Message	For this setting, press the measurement button to go to the print message input screen. You can enter print messages in this screen.	

Settings for Patient ID Input Interface

In the patient ID input interface menu, you can perform settings related to the position of the start byte for reading data, number of bytes for reading input data, method for inputting data, and connection tests for a USB device for inputting the patient ID.

How to Enter the Patient ID Input Interface Menu

MENU screen >  Input tab > [Input]

Setting items	Description	Factory settings
Update	Set the input methods to handle a USB device for automatic input (magnetic card reader, barcode reader, etc.) and a USB device for manual input (numeric keypad, etc.). [OVERWRITE]: Input data overwrites existing data. Select when a USB device for automatic input is connected. [ADD]: Input data is added to existing data. Select when a USB device for manual input is connected.	ADD
Start Position	Set the position of the byte at which to start reading data that is input from a USB device. You can only do this setting if [OVERWRITE] is selected for [Update]. [0 to 255]: Position of byte at which to start reading	0
Data Length	Set the number of bytes to read of the data that is input from a USB device. The length of data (number of bytes) read from the reading start byte position set for [Start Position]. You can only do this setting if [OVERWRITE] is selected for [Update]. [1 to 64]: Length of data to read	1
Connection Check	For this setting, press the measurement button to do a connection test for the USB device and display the results. If the connection is normal, a message to check at USB connection check screen appears. If the connection is not normal, an error message appears. [YES]: Go to the USB connection check screen. [NO]: End the USB connection test.	

Settings for Interface of Output of Measurements

You can send measurement results to an external device via LAN or RS-232C.

CAUTION

- Always turn off the power before connecting cables.
- Use the computer and other equipment that conform to the system standard IEC 60601-1 or IEC 60950-1 for the Full Auto Tonometer TX-20. Be sure that the entire system conforms to IEC 60601-1-1. Be sure to also use an isolation transformer conforming to IEC 60601-1 when a computer conforming to IEC 60950-1 is used.
Otherwise, electric shock may result. Ask your sales representative or local Canon dealer to connect it to a computer.
- Be sure to attach the provided ferrite core before connecting the LAN cable to the LAN connector.



Using the RS-232C or LAN connector

Ask your sales representative or local Canon dealer to check the connectable computer, connecting cables, setting and operation methods of the TX-20.

How to Enter the Output Settings Menu

MENU screen > Output tab > [Output]

Setting items	Description	Factory settings
Interface	Set the interface for external output. [OFF]: Turns off external output. [LAN]: Uses LAN interface. [RS-232C]: Uses RS-232C interface.	OFF

Settings for RS-232C

You can set the data format, transmission speed, transmission bits, parity, and stop bit in the RS-232C menu.

How to Enter the RS-232C Menu

MENU screen > Output tab > [RS-232C]

Setting items	Description	Factory settings
Format	Select the format to output measurement data. [T-2(EACH)]: T-2 format (1 data) [T-2(ALL)]: T-2 format (all data) [TX-F]: TX-F format [JMOIA]: JMOIA format [TX-20]: TX-20 format	TX-20
Baud rate	Select transmission speed. [300/600/1200/2400/4800/9600/19200]: Transmission speeds	9600
Bit	Select transmission bit. [7/8]: Transmission bit	8
Parity	Select parity. [ODD/EVEN/NONE]: Parity	NONE
Stop Bit	Select stop bit. [1/2]: Stop bit	1

Settings for LAN

You can set the DHCP, IP address, subnet mask, default gateway, user name, password, and IP address of the connected computer, DNS server, and name of data folder. You can also enable transmission of image data and do connection tests.

How to Enter the LAN MenuMENU screen >  Output tab > [LAN]

Setting items	Description	Factory settings
Stylesheet Output	Set to output stylesheets to the connected PC. [ON]: Output. [OFF]: Not output.	ON
DHCP	Set to enable DHCP. [ON]: Enable [OFF]: Do not enable	OFF
IP Address	Enter the IP address. [xxx.xxx.xxx.xxx]: IPv4 address	192.168.000.001
Subnet Mask	Enter the subnet mask. [xxx.xxx.xxx.xxx]: Subnet mask	255.255.255.000
Default Gateway	Enter the default gateway. [xxx.xxx.xxx.xxx]: IPv4 address	000.000.000.000
DNS Server	Input the DNS server. [xxx.xxx.xxx.xxx]: IPv4 address	192.168.000.100
User	Enter the user name of the computer that is connected. For this setting, press the measurement button to go to the user name input screen. Supported character code: ASCII Maximum number of characters: 16	ftp
Password	Enter the user password of the computer that is connected. For this setting, press the measurement button to go to the password input screen. If the password is input, [*****] appears. Supported character code: ASCII Maximum number of characters: 16	ftp
Computer Name	Enter the name or IP address of the computer that is connected. For this setting, press the measurement button to go to the computer name input screen. Supported character code: ASCII Maximum number of characters: 64	192.168.000.010
Output Folder	Enter the path to a shared folder on the computer that is connected for saving measurements. Do not input "/" to indicate the root directory. For this setting, press the measurement button to go to the output folder input screen. Supported character code: ASCII Maximum number of characters: 31	Empty
Connection Check	For this setting, press the measurement button to do a connection test and display the results.	

Settings for Patient No.

You can set patient No., patient No. display, patient No. count up, output at measurement error, reset, in the Patient No. menu screen.

The patient No. is a number to identify patients. You can set any number from 000001 to 999999 and you can set the number to increase by one each time a measurement is done. Patient No. can be checked on the **Measurement** screen, **MENU** screen, printed data, and data output to external devices.

How to Enter the Patient No. Menu

MENU screen >  Patient tab > [Patient No.]

Setting items	Description	Factory settings
Display	Set the display of the patient No. on the Measurement screen and whether to increase the numbers in series. Increase the number after measurement when the PRINT button or the CLS key is pressed. [ON]: Displays. Numbers increase serially. [OFF]: Does not display. Numbers do not increase serially.	ON
Count up on Error	Set whether to increase the patient number serially even if the measurement results are just an error. [ON]: The patient No. increases serially. [OFF]: The patient No. does not increase serially.	ON
No.	Set the patient No. [000001 to 999999]: Six-digit number	1
Reset No.	Set the patient No. to "000001". For this setting, press the measurement button to display the message screen and you can select reset. [OK]: Resets. [CANCEL]: Does not reset.	

Settings for Patient ID

You can set the input screen display, number of characters to input, requirement to input, and type of characters to input, in the patient ID menu.

The patient ID is a code to identify patients. You can set up to 64 alphanumeric characters and symbols. The set patient ID can be checked on the **Measurement** screen, **PATIENT ID** screen, printed data, and data output to external devices.

How to Enter the Patient ID Menu

MENU screen >  Patient tab > [Patient ID]

Setting items	Description	Factory settings
Redisplay In Scn	Set whether to display the PATIENT ID screen after printing (PRINT button is pressed) or clearing data (CLS key is pressed). [ON]: Displays. [OFF]: Does not display.	OFF
Limit Data Length	Set whether to limit the number of characters in a patient ID. [OFF]: No limit. Maximum number of characters is 64. [1 to 63]: Limited. It is not possible to input more characters than the number set.	OFF
Always Input	Set whether to require that the patient ID be input. [ON]: Input is required. A message appears after the PATIENT ID screen is closed, and before the first printout, and before transmitting data. [OFF]: Input is not required.	OFF
Numeric	Set whether to display numbers in the character list on the PATIENT ID screen. [ON]: Displays. [OFF]: Does not display.	ON
Alphabet (Upper)	Set whether to display upper case letters of the alphabet in the character list on the PATIENT ID screen. [ON]: Displays. [OFF]: Does not display.	ON
Alphabet (Lower)	Set whether to display lower case letters of the alphabet in the character list on the PATIENT ID screen. [ON]: Displays. [OFF]: Does not display.	ON
Symbol	Set whether to display symbols in the character list on the PATIENT ID screen. [ON]: Displays. [OFF]: Does not display.	ON

Setting the Time, Date, and Power Saving Mode

You can set the date, time, and power saving mode, in the Other Menu.

On the TX-20, you can set the data and time so that the dates of measurements are output as data. Also, you can set the monitor to turn off and enter power saving mode if no operations are done for a set period.

Press a button on the operation panel or a soft key to recover from the power saving mode.

How to Enter the Other Menu

MENU screen >  Other tab > [Date & Time]

Setting items	Description	Factory settings
Date	Set the date. [2010 to 2099]: Year [1 to 12] or [JAN to DEC]: Month. The [Format] setting is [MDY]; if [DMY] is selected, then [JAN to DEC]. [1 to 31]: Date	
Time	Set the time. [00 to 23]: Hour [00 to 59]: Minute	
Format	Select the format to display the date. [YMD]: Year, month, date [MDY]: Month, date, year [DMY]: Date, month, year	MDY
Sleep Time	Set the time before entering power saving mode. [OFF]: Does not enter power saving mode. [5min]: Five minutes [10min]: Ten minutes [15min]: Fifteen minutes	5min

Settings for the LCD

You can set the brightness of the LCD screen in the LCD menu.

How to Enter the LCD Menu

MENU screen >  Other tab > [LCD]

Setting items	Description	Factory settings
Adjust Brightness	Adjusts the brightness of the LCD. For this setting, press the measurement button to go to the LCD brightness adjustment screen. [1 to 100]: LCD brightness	50

Checking Firmware

You can check the name and version of the firmware in the firmware menu.

How to Enter the Firmware Menu

MENU screen >  Other tab > [Firmware]

Setting items	Description	Factory settings
Name	Displays the name of the firmware.	TX-20
Version	Displays the version of the firmware.	

6 Maintenance

This section outlines the maintenance instructions for the TX-20: daily inspection, cleaning, disinfecting, and refilling chin rest paper.

CAUTION

- For safety reason, before using the TX-20, be sure to perform the daily inspection.
- Have a periodically inspection performed for the TX-20 at least once a year by a Canon designated representative to maintain its performance and reliability.

Daily Inspections

Perform the following inspections before using the TX-20 to ensure that it is used safely and correctly. If a problem is found during the inspection and you are unable to correct the problem, please contact your sales representative or local Canon dealer.

Checks before turning on the power

Check the following items before turning on the power.

Cables

- 1) The power cord and connection cable are not damaged and their insulation is not torn.
- 2) The power cord is fully and securely inserted into the AC connector on the main unit and the AC outlet.

Main unit

- 1) There are no scratches or dirt on the nozzle and objective lens. Clean the nozzle and objective lens if they are dirty (see page 72).
- 2) There is no dirt on the monitor. Clean the monitor if it is dirty.
- 3) The forehead rest is disinfected (see page 73).
- 4) The chin rest paper is loaded (see page 74).
If the chin rest paper is not being used, disinfect the chin rest (see page 73).
- 5) The printing paper is loaded (see page 74).
- 6) The covers and other parts are neither damaged nor loose.

Checks after turning on the power

Check the following items after turning on the power.

- 1) The power lamp lights.
- 2) The chin rest moves up and down smoothly as the CHIN REST button is pressed.
- 3) The measurement unit moves smoothly in all directions (forward, back, left, right, up, and down) as the operation lever and the measurement unit vertical movement ring are operated.

Cleaning and Disinfection

CAUTION

Always disinfect the nozzle after measuring an infected patient's eyes

A patient's tears may get on the nozzle and cause a secondary infection. Always use a cotton swab or other appropriate material that has been treated with disinfectant to wipe the nozzle.

Nozzle

When dust or dirt is stuck on the tip of or inside the nozzle, press the CHECK button to blow it off.

Optical Component of Nozzle




- Do not wipe or rub the optical component if there is dirt or dust on it. Doing so could scratch the surface of the optical component.
- Do not wipe the optical component with disinfectant ethanol, eyeglass cleaner, or silicone-coated paper. Doing so could corrode the surface of the optical component or leave streaks. Also, wiping with too much lens cleaner causes streaking.
- For information about lens cleaning paper and lens cleaner, contact your sales representative or local Canon dealer.
- Be sure to put the cap over the optical component of the nozzle when the TX-20 is not in use.

Clean the optical component of the nozzle as indicated below when dust is on it, or if it is stained with tears.


1 Blow away any dust or dirt.

Remove the brush from the blower brush and blow off the dust on the optical component.

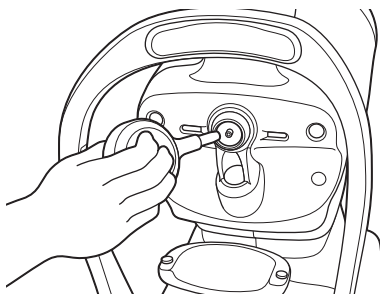

2 Wipe the optical component.

Slightly dampen a lens cleaning paper with lens cleaner and lightly wipe the optical component. Wipe it again if it needs it.

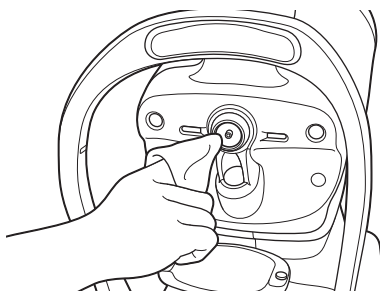



Use lens cleaning paper and lens cleaner designated by Canon.

Forehead Rest

CAUTION

- To prevent the risk of infection, wipe the forehead rest with disinfectant ethanol for each patient.
- The forehead rest may be corroded if a disinfectant other than ethanol is used.
- For details on how to disinfect, consult a specialist.

Use a sanitized gauze or wipe that includes disinfectant to clean the forehead rest for each patient.

Exterior and monitor screen of the TX-20

WARNING

- Before cleaning the exterior of the TX-20, be sure to turn off the power of all the devices, and unplug the power cord from the AC outlet.
- Do not use any flammable solvent such as alcohol that is not for cleaning, benzene, or thinner. Otherwise, fire or electric shock may result.

Clean the exterior, rollers, and monitor screen as follows if they are dirty.

1 Turn off the power.

Turn off the power switch and disconnect the power cable from AC receptacle.

2 Wipe the exterior and rollers.

Wipe them with a soft cloth that has been soaked in diluted neutral detergent and firmly wrung out. Then, wipe the exterior and rollers again with a cloth that has been soaked in water and firmly wrung out.

3 Wipe the monitor screen.

Wipe it lightly with a soft cloth such as an absorbent cotton rag dampened with a small amount of ethanol.

Refilling Chin Rest Paper

CAUTION

If the chin rest paper is not being used, be sure to disinfect the chin rest whenever the patient is changed. For details on how to disinfect, see the Forehead Rest section.

- 1 Pull out the right and left chin rest holding pins.
- 2 Insert the holding pins into the holes on the right and left of the chin rest paper.
- 3 Attach the chin rest paper to the chin rest with the pins pressing into the holes on the chin rest.



Insert the pins straight into the holes.



Chin rest paper is a consumable product (sold separately). To purchase chin rest paper, contact your sales representative or local Canon dealer.


Replacing Printing Paper

Replace the roll of printing paper as soon as possible after the red line appears on the paper.



Do not replace the paper unless the TX-20 is correctly installed in a flat location as specified. Changing the paper while the TX-20 is not level may cause a paper jam.

- 1 Open the paper cover.
Pressing the paper cover button releases the lock and raises the paper cover.

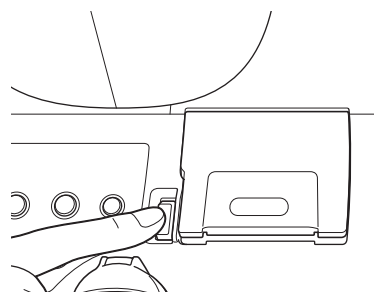

CAUTION

Do not touch the area around the printer or open the paper cover during or immediately after printing

The printer's thermal head and the surrounding area are very hot during and immediately after printing. To prevent burns, wait for the printer to cool before opening the paper cover to replace paper.

2 Remove the remaining paper.

Remove the old roll of paper.

**CAUTION****Do not touch the cutter for the printer**


Do not touch the cutter for the printer. Similarly, instruct the patient not to touch the cutter for the printer.

Otherwise, it may result in injury.

3 Put new paper into the printer.

Load the roll of paper into the printer as illustrated.


Be sure to check that the direction of roll is the same as shown in the illustration.

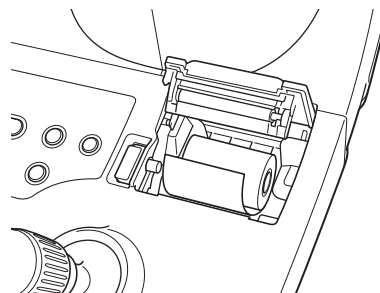**CAUTION**


Print quality, thermal head service life, cutting operations, and auto cutter service life may not be guaranteed if thermal paper not specified by Canon is used. Always use thermal paper specified by Canon.

4 Set the printing paper.

Pull out about 5 cm of the leading edge of the paper through the paper slot while the paper cover is open. Set the printing paper straight so it does not skew.

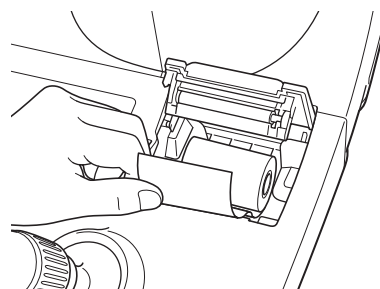**5** Close the paper cover.

When the paper is set correctly, close the paper cover until it is in the lock position. The paper automatically feeds and is cut.

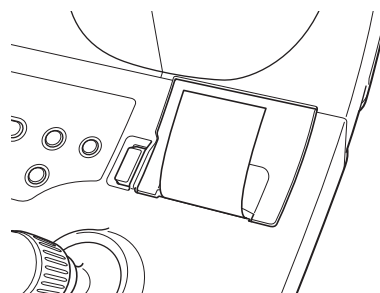**CAUTION**

After closing the paper cover, do not pull the thermal paper out before the paper is cut automatically. There is a risk of damaging the printer.



If the paper is fed through the printer incorrectly, it may cause a paper jam. If the paper jams, open the paper cover and adjust the paper so it feeds correctly.

6 Cut the paper.

Paper that is automatically cut has 1.5 mm of uncut paper in the middle of the cut so it does not fall off. Hold the tip of the paper and slowly pull it to cut it.



Always remove the automatically cut paper.

Printing again while the cut paper is still attached may cause paper jams.

Preparations for Packing and Transport

Before packing and transporting the TX-20, move its measurement unit and chin rest to the optimal position for packing and moving.

The measurement unit moves to the central lowest position. The chin rest moves all the way down.

1 While pressing the R/L button, turn on the power.

Continue pressing the R/L button after turning on the power until the buzzer sounds.

The chin rest and the measurement unit automatically move to the packing and transport position.

The following message appears on the monitor.

Transportation Package Position.
Please power down and remove the printing paper.

After that, the TX-20 cannot be operated. To do measurements again, the power must be turned off and on again.

2 Turn off the power.

Remove the paper from the printer after turning off the power. Unplug the power cord from the electric outlet and remove the cables to any connected devices.

7 Troubleshooting

If any of the symptoms listed below occur or if a message appears while using the TX-20, try the relevant remedy described below.

If these remedies do not eliminate the symptoms or the message, turn off the power and contact your sales representative or local Canon dealer. Be ready to provide a detailed description of the symptoms and message.


Symptom and Remedy

Symptom	Cause	Remedy
Measurements are not stable.	Eyelashes or eyelid are blocking the pupil.	Tell the patient to open his/her eye wide. Or, carefully lift the upper eyelid and perform the measurement.
	The patient's eye is not still.	Relax the patient and have him/her look at the eye fixation lamp.
	External light is getting in.	Do not allow bright light or sunlight to shine directly on the patient's eye.
Too few lines for the received light level mark.	Eyelashes or eyelid are blocking the pupil.	Tell the patient to open his/her eye wide. Or, carefully lift the upper eyelid and perform the measurement.
	Blinking is increasing.	Tell the patient to stop blinking or carefully lift the upper eyelid and perform the measurement.
	Patient's eye has low corneal reflectivity.	In this case, the received light level marks do not increase.
Nothing appears on screen even though the power is on and the POWER lamp is flashing.	Power saving mode.	Press a button on the operation panel, a soft key on the monitor panel, or the measurement button, or tilt the operation lever.
Paper jam occurs when printing.	<ul style="list-style-type: none"> The printing paper is fed through the printer incorrectly. The printing paper is set incorrectly. The printing paper from last use has not been removed. 	<ol style="list-style-type: none"> 1) Remove printing paper if there is any left from the last use. 2) Open the paper cover. 3) Check if the printing paper is fed through the printer aslant. 4) Take out the printing paper and remove the part of paper which is wrinkled or creased. 5) Set the printing paper correctly. 6) Check if printing is done correctly.
The paper cutter stops midway in a cut.	<ul style="list-style-type: none"> The printing paper is set incorrectly. Serious paper jam occurred. Abnormality occurred in printer. 	<ol style="list-style-type: none"> 1) If the message is displayed in the screen, record the message. 2) Press the paper cover button to move the cutter to its initial position. If pressing the button once does not move the cutter, press it several times. 3) Open the paper cover. 4) Check if the printing paper is fed through the printer incorrectly. 5) Take out the printing paper and remove the part of paper which is wrinkled or creased. 6) Set the printing paper correctly. 7) Check if printing and cutting are done correctly. 8) If this symptom occurs again, turn the power off and then on and check if printing is done correctly.

Message and Remedy

The message screen appears if a problem occurs or if there is a notification.
If this happens, check the displayed message, refer to the message codes and error codes on the following list and follow the instructions in the remedy. If the message screen does not close or if a message not shown in the table is displayed, contact your sales representative or local Canon dealer.

Message Screen

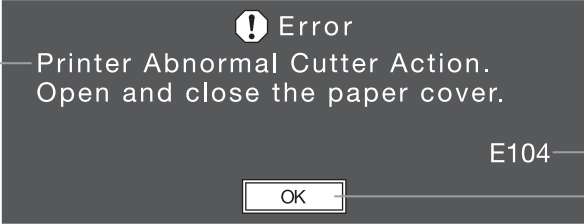

1 — Message content
Displays the content of the notification and the status of the problem that has occurred. It may also display the error code which indicates the details of the error status.

2 — Message code
The message category is indicated by the first letter.
Fxxx: The customer cannot resolve this error.
Exxx: The customer may be able to apply a remedy to resolve this error.
Wxxx: Alerts the customer to be careful.

3 — Check icon

Message Screen Operations

- 1 Check the message.

2 Select the check icon.
Rotate the measurement unit vertical movement ring to select the confirm icon, and press the measurement button.

3 Apply the remedy.
Apply the remedies as shown in the message.

Fatal

Code	Cause	Remedy
1	Abnormality occurred in transmission between boards.	1) Record the message code and contents (including error codes). 2) Turn the power off and then on again. If it does not return to normal, stop using the unit and contact your sales representative or local Canon dealer.
2	Abnormality occurred when CCD turned on.	
3	Abnormality occurred in light sensor for intraocular pressure measurement.	
4	Abnormality occurred in pressure sensor for intraocular pressure measurement.	
5	Abnormality occurred while checking memory.	
6	Abnormality occurred while writing to memory.	
7	Abnormality exists in data being stored.	
8	Abnormality occurred while measurement unit was detecting position during move to initial position.	
9		
10		
11		
12		
13		
14		
15		
16		
17		
18		
19		
20		
21		
22	Operation stopped because a button or the operation lever was operated while measurement unit was moving to initial position.	
23	Abnormality found on the charging function for air blower during the initial checking after power on.	

7 Troubleshooting

Error

Code	Cause	Remedy
100	Found abnormality in setting during initialization operation when power was turned on. Return to factory settings.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Turn the power off and then on again.
101	Abnormality occurred while checking clock.	
102	Failed to read time.	
103	Failed to set time.	
104	Abnormality occurred during printer cutter operation.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Open and close the paper cover and check if the error recurs. 3) If the error occurs again, turn the power off and then on.
105	The voltage being supplied to the printer is low.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Turn the power off and then on again.
106	Abnormality occurred at printer cutter position.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Open and close the paper cover and check if the error recurs. 3) If the error occurs again, turn the power off and then on.
107	Printer operation was not done normally.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Turn the power off and then on again.
108	Printer power was not turned on or off normally.	
119	Abnormality occurred while charging driver for air blower when CHECK button was pressed.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Press the CHECK button again and check if the error recurs.
120	Abnormality occurred while measurement unit was detecting position during move to initial position.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Turn the power off and then on again.
123	Abnormality occurred in USB transmission.	
124		
125	Abnormality occurred in LAN transmission.	
126		
127		
128		
129		
130		
131		
132		
133		
134		
135		
136		
137		
138		
139		
140		
141		
142		
143		
144		

Warning

Code	Cause	Remedy
201	The printer is overheating.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Wait a while and then print again.
202	The paper cover is open.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Close the paper cover.
203	There is no paper.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Replace the paper.
204	Printing stopped because a button was operated during printing.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Print again.
206	Operation stopped because a button was operated during printer initialization.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) No specific remedy is needed, you can do the next operations.
209	Operation stopped because a button was operated while sending data via RS-232C.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Do the transmission again.
211	The RS-232C connection is abnormal.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Check that the RS-232C cable is connected correctly (use a cross-cable for the RS-232C cable). Also, check the RS-232C transmission settings on the connected computer.
212	Abnormality occurred during RS-232C data transmission.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Check the RS-232C transmission settings.
213		
214		
219		
220		Check the RS-232C transmission settings in the software on the connected computer.
223	No USB device is connected.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Check that the USB device is connected correctly.
224	An unsupported USB device is connected.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Connect an HID class USB device. For details, contact your sales representative or local Canon dealer.
225		1) Record the message code, contents (including error codes) and then execute [OK]. 2) Connect a USB device that is compatible with HID class version 1.1. For details, contact your sales representative or local Canon dealer.
226		1) Record the message code, contents (including error codes) and then execute [OK]. 2) Connect a USB device that has a keyboard protocol code. For details, contact your sales representative or local Canon dealer.
227		1) Record the message code, contents (including error codes) and then execute [OK]. 2) Connect a USB device that has a low transmission speed. For details, contact your sales representative or local Canon dealer.
228	Two or more USB devices are connected.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Reduce the number of connected USB devices to one. For details, contact your sales representative or local Canon dealer.

7 Troubleshooting

Warning

Code	Cause	Remedy
229	Detected over current on the USB device connection.	<ol style="list-style-type: none"> 1) Record the message code, contents (including error codes) and then execute [OK]. 2) Check the connected USB device — it may not be compatible with the USB standard. Or, connect a USB device that is compatible with the standard. For details, contact your sales representative or local Canon dealer.
230	Could not do auto-alignment because could not find eye within set time period. The following problems may have occurred:	
	The patient's eye is outside the possible range for auto-alignment.	<ol style="list-style-type: none"> 1) Adjust the measurement unit and the chin rest to position the entire pupil of the patient's eye in the monitor. 2) Perform the measurement again.
	The patient's eye moved. Or the eye is nystagmic.	<ol style="list-style-type: none"> 1) Instruct the patient to stare at the eye fixation lamp, confirm the patient's eye has stopped moving and then perform the measurement. 2) If a nystagmus or another problem causes the patient's eye to move constantly during auto-alignment, press the measurement button when the corneal light is clearly visible. You can perform measurements during auto-alignment. However, the positioning alignment may not be accurate. Note that the measurement may not be accurate when measuring this way.
	Eyelid or eyelashes are blocking the pupil.	<ol style="list-style-type: none"> 1) If the eyelid or eyelashes are blocking the pupil, instruct the patient to keep his/her eyes wide open until the measurement is finished. 2) If the patient cannot open his/her eyes, the examiner should lightly hold the upper eyelid and do the measurement while helping to lift the patient's eyelid. Be careful to not hurt the patient.
	Mascara is blocking the pupil.	<ol style="list-style-type: none"> 1) Adjust so the corneal light point appears. 2) Do the measurement again.
	The pupil may be dystrophic or deformed. Or, there is a possibility the patient has undergone iris surgery.	<ol style="list-style-type: none"> 1) Adjust the measurement unit and the chin rest to position the corneal light point so it appears near the center. 2) Do the measurement again.
	A corneal disorder or some other abnormality exists on the surface of the cornea.	<ol style="list-style-type: none"> 1) Do the measurement in manual alignment mode. However, an error may occur even in manual alignment mode depending on the severity.
	External light is shining on the patient's eye.	<ol style="list-style-type: none"> 1) Do not allow bright light or sunlight to shine on the patient's eye during measurement.

Warning

Code	Cause	Remedy
231	Auto-alignment stopped because the limit of motion (to right) of the measurement unit was detected during auto-alignment.	<ol style="list-style-type: none"> 1) Record the message code, contents (including error codes) and then execute [OK]. 2) Be sure that the patient's forehead and chin are firmly against the forehead rest and chin rest and then do measurement again.
232	Auto-alignment stopped because the limit of motion (to left) of the measurement unit was detected during auto-alignment.	
233	Auto-alignment stopped because the limit of motion (up) of the measurement unit was detected during auto-alignment.	
234	Auto-alignment stopped because the limit of motion (down) of the measurement unit was detected during auto-alignment.	
235	Auto-alignment stopped because the limit of motion (forward) of the measurement unit was detected during auto-alignment.	
236	Auto-alignment stopped because the limit of motion (backward) of the measurement unit was detected during auto-alignment.	
237	Auto-alignment stopped because the set stop position of the nozzle was detected during auto-alignment.	<ol style="list-style-type: none"> 1) Record the message code, contents (including error codes) and then execute [OK]. 2) Do the following and do the measurement again. <ul style="list-style-type: none"> • Press the LIMITER button to reset the position where the nozzle stops. • Be sure that the patient's head and chin are firmly against the forehead rest.
238	Operation stopped because a button or the operation lever was operated during auto-alignment.	<ol style="list-style-type: none"> 1) Record the message code, contents (including error codes) and then execute [OK]. 2) Do the measurement again.
241	Measurement stopped because an intraocular pressure measurement error occurred during full auto-alignment measurement.	<ol style="list-style-type: none"> 1) Record the message code, contents (including error codes) and then execute [OK]. 2) Do the measurement again.
242	The measured intraocular pressure is over the value set for the [Upper IOP Warning].	
243	The measured intraocular pressure is below the value set for the [Lower IOP Warning].	
244	Auto-alignment stopped because the measurement unit and the chin rest may be interfering with each other.	<ol style="list-style-type: none"> 1) Record the message code, contents (including error codes) and then execute [OK]. 2) Be sure that the patient's head and chin are firmly against the forehead rest.

7 Troubleshooting

Warning

Code	Cause	Remedy
245	Intraocular pressure measurements are dispersed.	<ol style="list-style-type: none"> Record the message code, contents (including error codes) and then execute [OK]. Do the measurement again and check if there is a reason for the dispersion of the measured values. The following may be some reasons why the measured values are dispersed. <ul style="list-style-type: none"> Eyelashes or eyebrow are blocking the pupil. The patient's eye moved. Interference from external light.
248	Operation stopped because a button or the operation lever was operated while the measurement unit was moving left or right.	<ol style="list-style-type: none"> Record the message code, contents (including error codes) and then execute [OK]. Do the operation again.
249	Operation stopped because limit of motion was detected while the measurement unit was moving left or right.	<ol style="list-style-type: none"> Record the message code, contents (including error codes) and then execute [OK]. Move the measurement unit so the patient's eye appears in the monitor.
250	Operation stopped because a button or the operation lever was operated while the measurement unit was moving to the initial position.	<ol style="list-style-type: none"> Record the message code, contents (including error codes) and then execute [OK]. Do the operation again.
251	Patient ID has not been input.	<ol style="list-style-type: none"> Select [YES] to input the patient ID and then input the patient ID. Select [NO] if you are not going to input the patient ID. If the patient ID does not normally need to be input, on the patient ID menu set the [Always Input] setting to [OFF].
252	Confirmation about whether to delete all measurement data.	<ol style="list-style-type: none"> Select [OK] to delete all the measurement data. Select [CANCEL] to not delete all the measurement data.
253	Appears if you try to change the IP address while [DHCP] in the LAN menu was [ON].	<ol style="list-style-type: none"> Set the [DHCP] setting to [OFF] when setting the IP address (see page 65).
254	Appears if you try to change one of the following settings while [Display] in the patient No. menu is [OFF]. <ul style="list-style-type: none"> [No.] on the patient No. menu [Count up on Error] on the patient No. menu [Reset No.] on the patient No. menu 	<ol style="list-style-type: none"> You can change the setting by setting [Display] to [ON] (see page 67).
255	Appears if you try to change one of the following settings while [Update] in the patient ID input interface menu is [ADD]. <ul style="list-style-type: none"> Position of byte at which to start reading Number of bytes in the data to read 	<ol style="list-style-type: none"> You can change the setting by setting [Update] to [OVERWRITE] (see page 64).
257	Appears if you try to input more than 64 characters while [Update] in the patient ID input interface menu is [ADD].	<ol style="list-style-type: none"> Do one of the following remedies. <ul style="list-style-type: none"> Use the BackSpace key on the USB device to reduce the number of characters that have been input. Press the CLS key to delete all the characters that have been input. Press the OK key to end the connection test.
258	Confirm whether to save settings when changing the LAN settings.	<ol style="list-style-type: none"> Select [YES] to save the settings. Select [NO] to not save the settings.
259	Confirmation about whether to clear all input data.	<ol style="list-style-type: none"> Select [OK] to clear. If you do not want to clear the study results, select [CANCEL].

Warning

Code	Cause	Remedy
260	Confirmation for executing reset of patient No.	1) Select [OK] to reset. If you do not want to reset the study results, select [CANCEL].
261	Appears if the format of the text input for the computer name is incorrect.	1) Input it again in the correct format (see page 65).
262	Appears if the format of the text input for the user name is incorrect.	
264	Appears if the format of the text input for the path to the folder to store data is incorrect.	
265	Confirm whether to save settings when changing the LAN settings and then doing a LAN connection test.	1) Select [OK] to save the settings. 2) Select [CANCEL] to not save the settings and not do the connection test.
269	Operation stopped because a button was operated during LAN transmission.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Do the operation again.
271	Abnormality occurred in LAN transmission.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Do the transmission again.
272		
273		
280		
281		
282		
283		
284		
285	Abnormality occurred in LAN transmission. The user name is not correct.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Check that the [User] in the LAN menu is set correctly (see page 65).
286	Abnormality occurred in LAN transmission. The password is not correct.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Check that the [Password] in the LAN menu is set correctly (see page 65).
297	Abnormality occurred in LAN transmission. The file you are trying to send already exists.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) If a file in the folder already exists with the same name, change the name of one of the files. Or, use a different folder. 3) Do the transmission again.
298	Abnormality occurred in LAN transmission.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Do the transmission again.
299	Abnormality occurred in LAN transmission. The output folder does not exist.	1) Record the message code, contents (including error codes) and then execute [OK]. 2) Check that the [Output Folder] in the LAN menu is set correctly (see page 65). Or, check that a folder exists where you are outputting the data.

7 Troubleshooting

Warning

Code	Cause	Remedy
320	Because of a dead battery (for backup), automatic updating of date stops.	<p>1) Measurements can be done. Select either one from the followings and use the TX-20 as it is.</p> <ul style="list-style-type: none">• Select [OK]. Reflect the clock on the measurement results. Set the correct clock again on the displayed clock input screen. However, the clock you have set again is valid only when the power is turned on.• Select [CANCEL]. The clock is not reflected correctly on the measurement results. The clock input screen is not displayed and the TX-20 is reset at the initial value 2010/1/1. <p>If this error occurs, the clock function of the main unit is not kept. Therefore, you need to set the clock each time you turn on the power. We recommend replacing the battery. For replacing the battery, ask your sales representative or local Canon dealer.</p>
329	The USB device of full-speed specifications was connected to this product.	<p>1) Record the message code, contents (including error codes) and then execute [OK].</p> <p>2) Connection check is necessary. Ask your sales representative or local Canon dealer to perform the connection check.</p>

Appendix

Specifications

Intraocular pressure measurement	
Measurement type	Air-puff type
Operation distance	11 mm
Range of intraocular pressure measurement	0 to 80 hPa (0 to 60 mmHg)
30 mode	0 to 40 hPa / 33 to 80 hPa changes automatically (0 to 30 mmHg / 25 to 60 mmHg changes automatically)
60 mode	33 to 80 hPa (25 to 60 mmHg)
Displayed unit of intraocular pressure	1 hPa / regular measurements: 0.1 hPa (1 mmHg / regular measurements: 0.1 mmHg) * 1.33 hPa = 1 mmHg
Observation range	13.2 x 10 mm or more
Movement range	
Measurement unit	Front/back: 40 mm Right/left: 90 mm Up/down: 30 mm
Chin rest	60 mm
Printer	Thermal printer (auto cutter, paper width 58 mm)
Interface	USB-A: 1 port RS-232C: 1 port, LAN: 1 port (only one selectable)
Display	14.5 cm (5.7 inch) TFT color LCD / Tilttable
Power saving function	5 min. / 10 min. / 15 min. / OFF (4 levels)
Power supply rating	AC 100 V to 240 V, 50/60 Hz, 0.8 A to 0.4 A
Dimensions	260 (W) x 490 (D) x 500 (H) mm
Weight	15 kg

Stylesheet File

Setting the stylesheet to be output when transmitting via a LAN allows you to view the measurements on a web browser such as Windows Internet Explorer.

The stylesheet file defines the layout of the display when measurement files are displayed on a web browser.



- Set whether to output the stylesheet file with the [Stylesheet Output] setting on the LAN menu (see page 65).
- Use a web browser that supports XML.



Ask your sales representative or local Canon dealer to connect the LAN cables and do the settings, and to explain the data transfer and operation methods.

How to View Measurements on a Web Browser

1 Open the folder in which to store data.

Open the folder that was set on the [Output Folder] in the LAN menu for outputting data (see page 65). A maximum of two types of files can be transmitted at one time. The files that can be transmitted are shown below.

- **Measurement files**

Includes patient information and intraocular pressure values.

File name

Patient information (patient ID etc.)_output date (yyyymmdd)_output time (hhmmss)_

CANON^TX-20.xml

Example: 999999_20110204_174235_CANON^TX-20.xml

- **Stylesheet file**

File name

TX_Stylesheet.xsl

2 Open the measurement file on a web browser.

The patient information, measurements, and other information are displayed on the web browser. An example of a display is shown below on the following.



Measurement files that are created when [Stylesheet Output] is set to [OFF] cannot be viewed on a web browser even if a style sheet file is output later.

TX-20 Measurement Data

Company Name	CANON INC.
Model Name	TX-20
Serial No.	999999
Firmware Version	1.0.0.0
XML Version	1.2
Date	2011/02/03
Time	19:36:28
Patient No.	123456
Patient ID	ABCDEFG

Right Eye				Left Eye			
TONO	IOP[mmHg]	IOP[hPa]	Low Confidence	TONO	IOP[mmHg]	IOP[hPa]	Low Confidence
1	15	20		1	15	20	*
2	14	19		2	16	21	*
3	15	20		3	10	13	
STD	14.5	19.3	-	4	9	12	
				5	10	13	
				6	22	29	*
				7	11	15	
				8	19	25	*
				9	10	13	
				10	10	13	*
				STD	10.3	13.7	-

EMC (Electromagnetic Compatibility)

The TX-20 is designed and tested to comply with IEC 60601-1-2 (EN 60601-1-2) which are applicable regulations regarding EMC for medical devices and need to be installed and put into service according to the EMC information stated as follows.

If this equipment causes harmful interference to other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving device.
- Increase the separation between the equipment.
- Connect this device into an outlet on a circuit different from that to which the other devices are connected.

If the problem cannot be solved with the above measures, stop using this equipment and consult your sales representative or local Canon dealer.

Precautions on EMC

1. Medical electrical equipment requires special precautions regarding EMC and must be installed and put into service according to the EMC information provided in the manual.
2. Portable and mobile RF communications equipment can affect medical electrical equipment.
3. Information regarding the cable affecting EMC is as follows:

To maintain the optimum EMC performance, use only the designated cables.

Name	Type	Length	Remarks
AC Power Cord	Non-Shielded	3.0 m	Supplied
RS-232C cable	Shielded	1.5 m	Not supplied
USB Ten Key Board	TK-TCM003SV	0.5 m	Not supplied
LAN cable	Non-Shielded	50 m	Not supplied

4. The use of accessories, transducers and cables other than those specified, with the exception of transducers and cables sold by Canon sales representative or local Canon dealer as replacement parts for internal components, may result in increased emissions or decreased immunity of the TX-20.
5. The TX-20 should not be used adjacent to or stacked with other equipment; if adjacent or stacked use is necessary, the TX-20 should be observed to verify normal operation in the configuration in which it will be used.

Guidance and Manufacturer's Declaration for Electromagnetic Emission

The TX-20 is intended for use in the electromagnetic environment specified below.

The customer or the user of the TX-20 should ensure that it is used in such an environment.

Emission Test	Compliance	Electromagnetic Environment – Guidance
RF emissions EN 55011 CISPR11	GROUP 1	The TX-20 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electromagnetic equipment.
RF emissions EN 55011 CISPR11	Class B	The TX-20 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions EN IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions EN IEC 61000-3-3*	Complies	

* Not applicable to regions where the rated voltage is less than 220V.

Appendix

Guidance and Manufacturer's Declaration for Electromagnetic Immunity

The TX-20 is intended for use in the electromagnetic environment specified below.

The customer or the user of the TX-20 should assure that it is used in such an environment.


Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Electrostatic discharge (ESD) EN IEC 61000-4-2	±(2, 4, 6) kV contact ±(2, 4, 8) kV air	±(2, 4, 6) kV contact ±(2, 4, 8) kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst EN IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge EN IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines EN IEC 61000-4-11	<5% U_T (>95% dip in U_T) for 0.5 cycles 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles <5% U_T (>95% dip in U_T) for 5 sec	<5% U_T (>95% dip in U_T) for 0.5 cycles 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles <5% U_T (>95% dip in U_T) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the TX-20 requires continued operation during power mains interruptions, it is recommended that the TX-20 be powered from an uninterruptible power supply.
Power frequency (50/60Hz) magnetic field EN IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: U_T is the AC mains voltage prior to application of the test level.

Guidance and Manufacturer's Declaration for Electromagnetic Immunity

The TX-20 is intended for use in the electromagnetic environment specified below.

The customer or the user of the TX-20 should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Conducted RF EN IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the TX-20, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter Recommended Separation Distances $d = 1.2 \sqrt{P}$ $d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey ^a , should be less than the compliance level in each frequency range ^b . Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF EN IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

- a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered.
If the measured field strength in the location in which the TX-20 is used exceeds the applicable RF compliance level above, the TX-20 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the TX-20.
- b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Recommended Separation Distances between Portable and Mobile RF Communications Equipment and Devices or the TX-20

The TX-20 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the TX-20 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the TX-20 as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Warranty and Repair Service

Service life

The service life of this product is eight years if specified inspections and maintenance are performed.

About repairs

If a problem cannot be solved even after taking the measures indicated in “7 Troubleshooting” (see page 77), contact your sales representative or local Canon dealer for repairs.

When requesting repair, please provide the following information.

Name of the instrument: TX-20

Serial number: 6-digit number on the rating label

Description of malfunction: Report as much detail as possible.

Time limit for supplying performance parts for repair

Performance parts (parts for repairs to maintain performance) will be stocked for eight years after production of the TX-20 is discontinued.

Expendable parts replaced by service personnel

The following parts cannot be replaced by the user. If these parts are found to be worn out or to have deteriorated during daily or regular inspection, contact your sales representative or local Canon dealer for repairs.

- Backup battery for the clock (type name: CR1632, battery life: approximately 8 years)







L-IE-5191C

Canon



Manufacturer:

CANON INC. Medical Equipment Group

30-2, Shimomaruko 3-chome, Ohta-ku, Tokyo 146-8501, Japan

Telephone: (81)-3-3758-2111

Distributed by:

CANON U.S.A., INC. Healthcare Solutions Division

15955 Alton Parkway, Irvine, CA 92618, U.S.A.

Telephone: (1)-949-753-4160



CANON EUROPA N.V. Medical Products Division

Bovenkerkerweg 59, 1185 XB Amstelveen, The Netherlands

Telephone: (31)-20-545-8926

CANON SINGAPORE PTE. LTD. Medical Equipment Products Division

1 HarbourFront Avenue, #04-01 Keppel Bay Tower, Singapore 098632

Telephone: (65)-6799-8888

CANON AUSTRALIA PTY. LTD. Optical Products Division

1 Thomas Holt Drive, North Ryde, Sydney N.S.W. 2113, Australia

Telephone: (61)-2-9805-2000

L-IE-5191C

0713P0.001

© CANON INC. 2011

Printed in Japan